

# CRISPR/Cas9 Genome Editing

## Application Guide

Introduction .....	5
I. Genome-wide CRISPR gene knockout kit.....	5
HDR-mediated CRISPR knockout kits.....	6
Package contents.....	6
Related Optional Reagents .....	7
Related OriGene Products .....	7
Notice to purchaser .....	7
Product Description.....	7
Experimental Protocol .....	9
KN2.0, non-homology mediated CRISPR knockout kits .....	13
Package contents.....	13
Related Optional Reagents .....	13
Related OriGene Products .....	13
Notice to purchaser .....	13
Product Description.....	13
Experimental Protocol .....	14
II. Precut All-in-one CRISPR vector cloning kits (pCas-Guide (SKU: GE100001), pLenti-Cas-Guide (SKU:GE100009) .....	18
Package contents .....	18
Related Optional Reagents.....	19
Related OriGene Products.....	19
Notice to purchaser.....	19
Production and Quality Assurance.....	19
Product Description .....	20
Experimental Protocols.....	20

I. Design target sequence .....	20
II. Addition of extra bases to the ends of the target sequence .....	20
III. Cloning the double-stranded oligos into the pCas-Guide vector .....	21
Lenti-based protocols:.....	22
III. Circular All-in-one CRISPR vectors (SKU: GE100002, GE100018, GE100022, GE100010, GE100045).....	24
Package contents .....	24
Related Optional Reagents.....	24
Related OriGene Products.....	24
Notice to purchaser.....	24
Product Description .....	25
Experimental protocol.....	26
IV. gRNA only Vectors, pGuide (SKU GE100042), pGuide-EF1a-GFP (SKU GE100044), pLenti-Guide-Puro (SKU GE100032).....	27
Package contents .....	27
Related Optional Reagents.....	27
Related OriGene Products.....	27
Product Description .....	27
Experimental protocol.....	28
V. Cas9 only vectors, .....	28
pLenti-Cas9 (SKU GE100028), .....	28
pLenti-Cas9-IRES-Puro (SKU GE100029),.....	28
pLenti-EF1a-Cas9-IRES-Puro (SKU GE100030), .....	28
pLenti-Cas9-P2A-tGFP (SKU GE100031).....	28
pAAVS1-Cas9-Puro-DNR (SKU GE100037).....	28
pAAVS1-Cas9-BSD-DNR (SKU GE100039).....	28
Package contents .....	28
Related Optional Reagents.....	28
Related OriGene Products.....	29

Product Description .....	29
Experimental Protocol for Lenti Cas9 Vectors .....	31
VI. T7 driven CRISPR/Cas system-pT7-Guide-IVT (SKU: GE100025) and pT7-Cas9 (SKU: GE100014) .....	31
Package contents .....	31
Related Optional Reagents.....	31
Related OriGene Products.....	32
Product Description .....	32
Experimental protocol .....	33
I. Design genomic target sequence and cloning into pT7-Guide-IVT vector .....	33
II. Producing gRNA and Cas9 mRNA using T7 in vitro transcription kits.....	34
VII. Cre expression vector for Cre-Lox recombination, pCMV6-Entry Cre (SKU: GE100017).....	34
Package contents .....	34
Related Optional Reagents.....	35
Related OriGene Products.....	35
Product Description .....	35
VIII. Cas9 D10A nickase vectors, pCas-Guide-Nickase (SKU GE100019), pT7-Cas9-Nickase (SKU GE100020).....	36
Package contents .....	36
Related Optional Reagents.....	36
Related OriGene Products.....	36
Product description .....	36
Experimental protocol .....	36
IX. CRISPR scramble controls, pCas-Scramble (SKU GE100003) and pCas-Scramble-EF1A-GFP (SKU GE100021).....	37
Package contents .....	37
Product description .....	37
X. Transgene Knockin via CRISPR at AAVS1 and ROSA26 Loci.....	38
Introduction.....	38


AAVS1 Transgene knockin via CRISPR.....	39
AAVS1 CRISPR/gRNA vector, pCas-Guide-AAVS1 (SKU GE100023) .....	39
AAVS1 donor vectors .....	41
AAVS1 Transgene knockin vector kit, puro (SKU GE100027, GE100047) or BSD (GE100036, GE100049) .....	45
pAAVS1-RFP-DNR (SKU GE100026), positive control.....	47
Cas9 insertion at AAVS1 locus .....	52
AAVS1 Cas9 insertion vector kits, Puro (SKU GE100038) and BSD (SKU GE100040).....	53
ROSA26 Transgene knockin via CRISPR .....	54
ROSA26 gRNA/CRISPR construct, pCas-Guide-ROSA26 (SKU GE100050) .....	54
ROSA26 donor vector, pROSA26-Puro-DNR (SKU GE100051),.....	55
ROSA26 Transgene knockin vector kit, puro (SKU GE100052) .....	57
FAQ.....	59

## Introduction

Cas9 based genome editing has become a popular tool for targeted genome manipulation because of its simplicity and high cutting efficiency. This system requires a functional cas9 protein and a guide RNA for effective double-stranded breakage at a desired site. OriGene has developed many CRISPR vectors, including All-in-one vectors which contain both guide RNA and Cas9 expression, T7 vectors and gRNA and Cas9 separate vectors. OriGene also designed a set of donor cassettes for HDR-based (Homology Directed Repair) donor vector construction, including Luciferase-Loxp-Puro-Loxp, GFP-Loxp-Puro-Loxp, RFP-Loxp-BSD-Loxp and BFP-Loxp-Neo-Loxp. Homologous arm sequences can be cloned flanking the donor cassettes.

OriGene CRISPR products also includes CRISPR gene knockout kits, gene harbor transgene insertion via CRISPR, synthetic gRNA.

Figure 1. Flow chart of CRISPR genome editing using HDR.

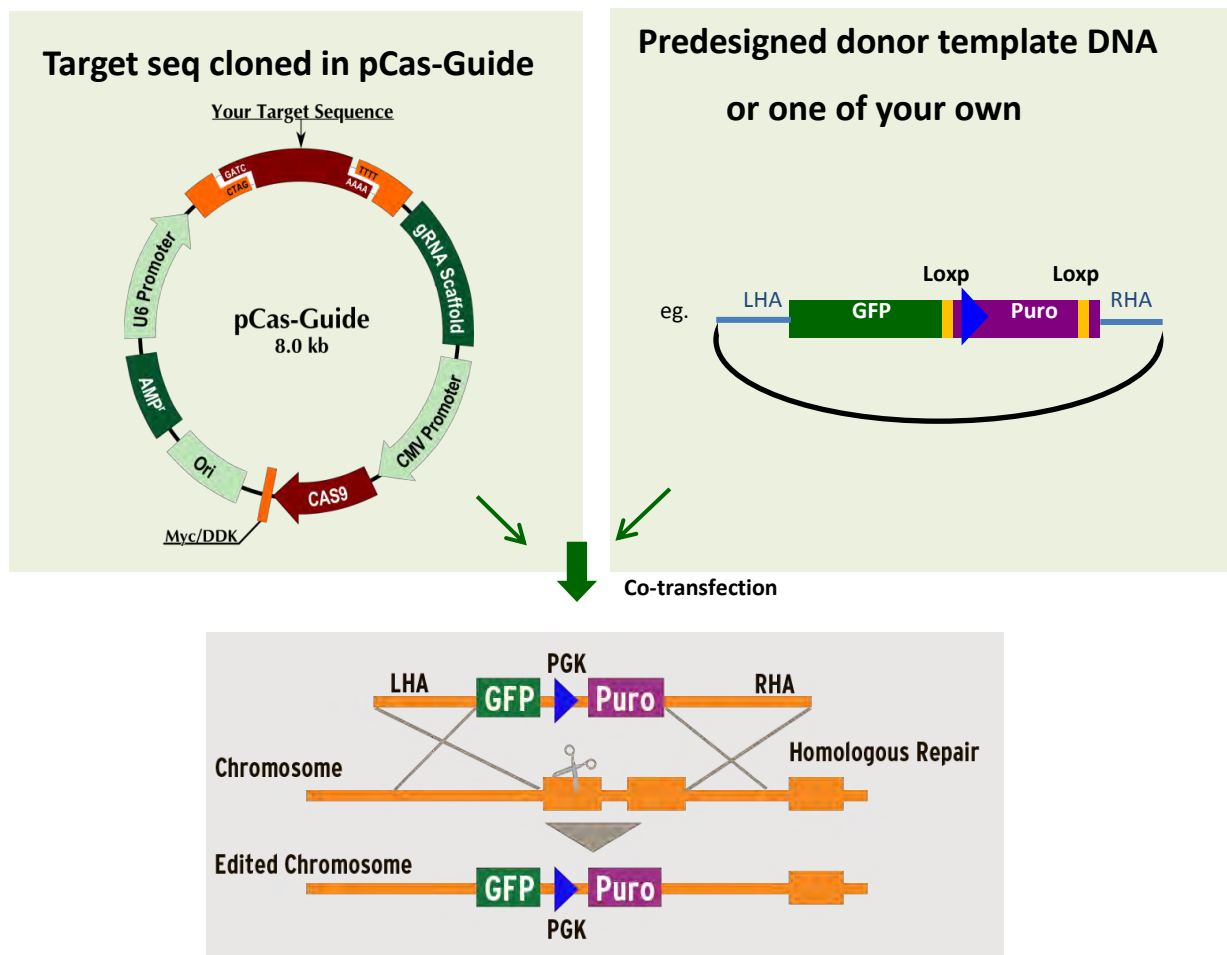

## I. Genome-wide CRISPR gene knockout kit

OriGene offers genome-wide gene knockout / knockin kits via CRISPR (human and mouse); one specific kit for each gene locus. Each kit contains 2 gRNA vectors and 1 donor DNA. The gene knockout/knockin kit is a complete kit to knockout any coding gene and knockin a selection cassette.

gRNA vectors are provided in pCas-Guide vector with a target sequence cloned. Both of the target sequences are located at the 5' end of the ORF; therefore gRNA vectors will make a precise cleavage at the 5' end of the ORF of the gene loci.

There are two types of CRISPR knockout kits: HDR mediated & KN2.0 non-homology mediated knockout kits. The difference of the two types kits is repair mechanism that is used to achieve gene knockout, Homology Directed Repair (HDR) and Non-Homology repair respectively.

**Table 2. Comparison of CRISPR gene knockout kits**

	HDR mediated	KN2.0 non-homology mediated
<b>gRNA vectors</b>	pCas-Guide	pCas-Guide
<b>Donor</b>	Donor cassette flanked by homologous arms	Donor cassette without homologous arms
<b>Cell spectrum</b>	Dividing cells	Dividing and non-dividing cells
<b>Knockout efficiency</b>	Medium	High

### ***HDR-mediated CRISPR knockout kits***

#### **Package contents**

- 2 vials of gRNA vectors, (SKU KNxxxxxxG1, KNxxxxxxG2), 3-5 µg DNA in TE buffer
- 1 vial of donor vector containing left and right homologous arms and a GFP-puro functional cassette (SKU KN2xxxxxD), 3-5 µg DNA in TE buffer
- 1 vial of negative scramble control vector (SKU GE100003), lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O, final concentration 100 ng/ µL.
- Certificate of Analysis
- Application Guide: <https://www.origene.com/support/product-support/product-manuals>

*Note: The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.*

### Related Optional Reagents

- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)

### Related OriGene Products

- Transfection reagent: <https://www.origene.com/products/others/transfection-reagents>
- Reagents and supplies for immunoblots: user preferred. OriGene has a selection of antibodies and detection reagents that are available at <https://www.origene.com/products/antibodies>
- DNA purification reagents <https://www.origene.com/products/others/nucleic-acid-purification-kits>
- qPCR reagents <https://www.origene.com/products/gene-expression/qpcr>
- CRISPR/Cas9 products <https://www.origene.com/products/gene-expression/crispr-cas9>

### Notice to purchaser

This product is for research use only. Use in and/or for diagnostics and therapeutics is strictly prohibited. By opening and using the product, the purchaser agrees to the following: The plasmids may not be distributed, resold, modified for resale or used to manufacture commercial products without prior written approval from OriGene Technologies, Inc. If you do not agree to the above conditions, please return the UNOPENED product to OriGene Technologies, Inc. within ten (10) days of receipt for a full refund.


### Product Description

The gene knockout/knockin kit is a complete kit to knockout any coding gene and knockin a functional cassette containing a reporter and selection marker. The donor plasmid contains left homologous arm and right homologous arm flanking the donor cassette, therefore, the donor cassette will be integrated into the genome via homology-directed repair (HDR) mechanism. The reporter, such as GFP, will be under the native promoter; the puromycin resistant gene is under PGK promoter. gRNA vectors are provided in pCas-Guide vector with a target sequence cloned. Both of the target sequences are located at the 5' end of the ORF; therefore gRNA vectors will make a precise cleavage at the 5' end of the ORF of the gene loci. A negative scramble gRNA control is also provided.

The Applications are:


1. You can knockin GFP reporter or other reporters for your promoter study.
2. Knock-out genes at the chromosomal level.

Fig. 2. Scheme of genome-editing knockout kit

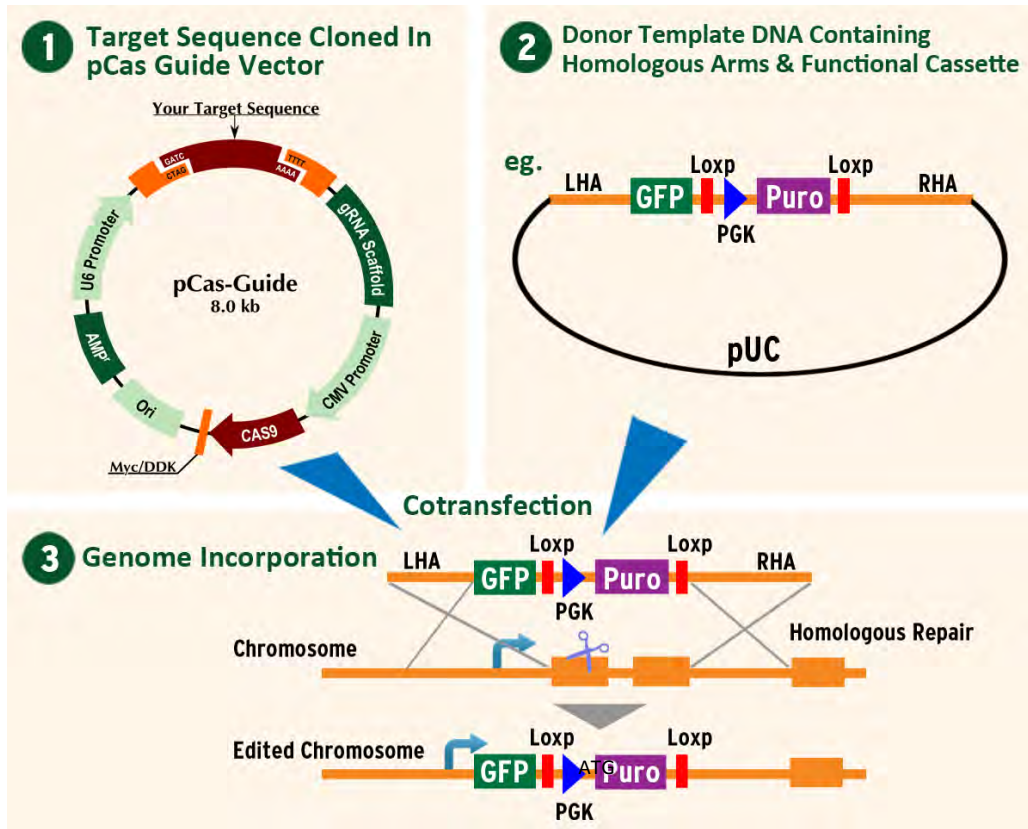

- 1 CRISPR/Cas cuts the double-stranded DNA at the targeting site
- 2 Donor template DNA provides the template for the homologous repair.
- 3 The functional cassette is incorporated into the genome when 1 + 2 are cotransfected.

Donor vector for each kit contains around 600 bp locus specific homologous sequence on each side of the donor selection cassette. LHA – left homologous arm, RHA – right homologous arm. Four different donor cassettes are offered for each gene:

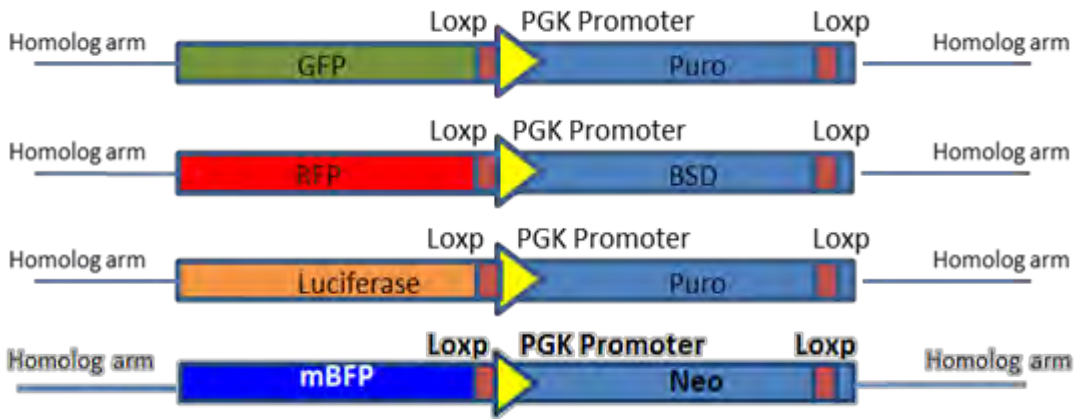




## Experimental Protocol

Each kit contains two gRNA vectors, one scramble negative control and one donor vector. To ensure high efficiency of cleavage, two gRNA constructs are provided. A scrambled control vector serves as the negative control.

A sample protocol listed below is for experiments performed in 6-well plates and using TurboFectin (cat# TF81001) as transfection reagent. If your experiments require a different size of culture plates, just scale up or down the reagents accordingly based on the relative surface area of your plate (Table 1). Different type of cells may need a different transfection reagent; please follow the manufacturer's corresponding protocol. OriGene just launched Virus-like Viromers which are best for difficult-to-transfect cells.

1. Approximately 18-24 hours before transfection, plate  $\sim 3 \times 10^5$  adherent cells in 2 ml culture medium into each well of a 6-well plate or  $\sim 5 \times 10^5$  suspension cells per well to obtain 50-70% confluence on the following day. The amount of cells varies depending on the size of your cells.
2. Transfection in a complete culture medium. Three separate transfections:

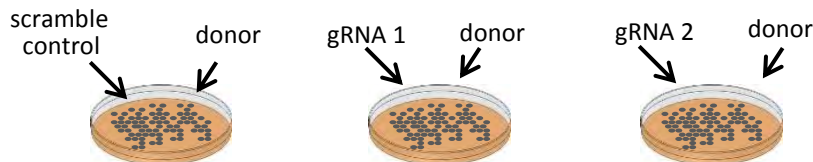

In a small sterile tube, combine the following reagents in the prescribed order. The order of reagent addition is important to achieve the optimal results.

- a. Dilute 1  $\mu\text{g}$  of one of the gRNA vectors (or scramble control) in 250  $\mu\text{L}$  of Opti-MEM I (Life Technologies), vortex gently. Then add 1  $\mu\text{g}$  of the donor DNA into the same 250  $\mu\text{L}$  of Opti-MEM I. Vortex gently. Two gRNA vectors and scramble control are in three separate tubes, so the gRNA efficiency can be tested individually.
- b. Add 6  $\mu\text{L}$  of Turbofectin 8.0 to the diluted DNA (not the reverse order) and pipette gently to mix completely.
- c. Incubate the mixture 15 minutes at room temperature.

**Note:** We recommend starting with the ratios of Turbofectin 8.0 and DNA listed in table 1; however, subsequent optimization may further increase the transfection efficiency.

- d. Add the mixture above drop-wise to the cells plated in step 1 (no need to change the media). Gently rock the plate back-and-forth and side-to-side to distribute the complex evenly.
  - e. Incubate the cells in a 5%  $\text{CO}_2$  incubator.
3. 48hr post transfection, split cells 1:10, grow additional 3 days; then split the cells again 1:10. Split cells 7 times in total. Since puromycin resistant gene in the donor vector

contains PGK promoter, the plasmid donor DNA before genomic integration will also provide puromycin resistance. The reason to grow cells for around 3 weeks before puromycin selection is to dilute out cells containing the donor as episomal form.

#### **Time lines of genome editing**

- ✓ CRISPR targeted gene knockout / knockin--- 1 week post transfection
- ✓ Episomal donor vector dilution with cell passaging--- 3 weeks post transfection

**Note 1.** *Since stable cell selection takes time, you can try to analyze the cells at P2 to detect genomic integration using genomic PCR (Fig 2). When designing primers for genomic PCR, one primer should be outside of the homologous arm region in donor DNA and one primer is in the functional cassette. Please see details in Fig. 3. The amplified PCR fragment is around 1kb. There could be some difficulties doing genomic PCR at this step before selection due to the percentage of edited cells and difficulties of genomic PCR. qPCR measuring the targeted mRNA level would not work due to the small percentage of edited cells.*

**Note 2.** *You might be able to use GFP to sort genomic edited cells between P2-P5 (Fig. 2). Since donor DNA contains 600bp left homologous arm sequence which is immediately upstream of ATG, Donor DNA transfected + scramble control could express weak or bright green fluorescence depending how much promoter sequence the left homologous arm sequence contains. The best case scenario is donor DNA + scramble gives weak GFP signal while after integration the promoter is strong and constitutive, so you can sort strong GFP positive cells, thus avoiding the lengthy donor vector dilution step before puro selection.*

4. Apply puromycin selection. Split P5 or P7 cells 1:10, then grow cells directly in the puromycin containing complete media in 10 cm dishes (apply puro selection at P6 might work, but you might get more random donor integrated clones than P8). The dose of puromycin needs to be determined with a kill curve to find out the lowest dose that kills the non-transfected cells completely 4-7 days post selection; the range of puromycin is 1µg/ml to 10µg/ml). Change the media every 2-3 days.

**Note:** *We recommend you still keep growing or freeze some of the transfected cells without selection; just in case, you need to perform the puromycin selection again.*

5. The puromycin resistant cells are ready to be analyzed for genome editing. WB can be used to measure gene knockdown if there is a good protein-specific antibody available (better after isolating individual cell colonies). You will also need to do genomic PCR to verify the integration of the functional cassette. You can directly sequence the amplified genomic fragment using the PCR primers to verify the integration in the genome.

**Note:** *Scramble control and donor vector will also give you some puromycin resistant cells as donor vector alone can randomly integrate into the genome too; however the efficiency should be a lot lower than with a specific gRNA. Therefore, you should get more colonies with gene specific gRNA than scramble control if the gene specific gRNA cleaves efficiently.*

6. Isolate individual cell colonies.

Two main methods, limiting dilution and cloning rings / cylinder.

1) Limiting dilution

This method is better to be used after puromycin selection. Dilute cells to seed about 1-2 cells/well in 96-well plate, after 1-2 weeks, observe under the microscope and select the wells only containing one cell colony, then further expand them to 6-well plate when they are confluent in the 96-well plate and so on.

2) Cloning rings / cylinder

This method can be used in the same time with puromycin selection. Seed cells at lower density, such as 5% confluency in a larger cell culture dish, such as 10cm dish, apply puromycin selection when seeding.

**Note 1: How to make bi-allelic knockout:** *If you isolate single cell colonies, in some cells gene knock-out may occur only in one allele; in some cells gene knock-out may occur in both alleles. If you only have monoallelic knockout (heterozygous) and you want to get bi-allelic knockout (homozygous), you can order another donor vector containing a different mammalian selection marker, such as blastocidin or neomycin resistant marker. Make sure the other allele is intact. You can confirm it with genomic PCR with a set of PCR primers amplifying the endogenous chromosome and sequence it. If this allele is targeted by Cas9/gRNA, repaired by NHEJ, introduced indels and the indels change the protein reading frame, then you have a bi-allelic knockout (one via HDR and one via NHEJ). If this allele is intact, you can do the knockout again. OriGene has both functional cassettes. You can do the knockout procedure again with the new donor vector to target the second allele (one allele is already targeted and replaced with GFP-puro cassette). Alternatively, you can use Cre (SKU GE100018) to flox out the puro cassette from your edited cells and use the same donor vector from the knockout kit to do the knockout again to target the second allele.*

**Note 2:** *If your gene is essential, you will not be able to get bi-allelic knockout. The solution is to do conditional knockout using LoxP system by introducing LoxP sites around the exon(s) to be knocked out.*

Table 1. Recommended starting transfection conditions for Turbofectin 8

Tissue Culture Vessel	Growth area, cm <sup>2</sup> /well	μg of DNA	Ratio of Turbofectin:DNA
96-well plate	0.35	0.1-0.15	3:1
24-well plate	2	0.5-1	3:1
12-well plate	4	1-2.5	3:1
6-well plate	9.5	1-5	3:1

35 mm plate	8	1-5	3:1
60 mm plate	20	2-10	3:1
100 mm plate	60	5-15	3:1

Figure 3. Diagram of cell passaging after transfection


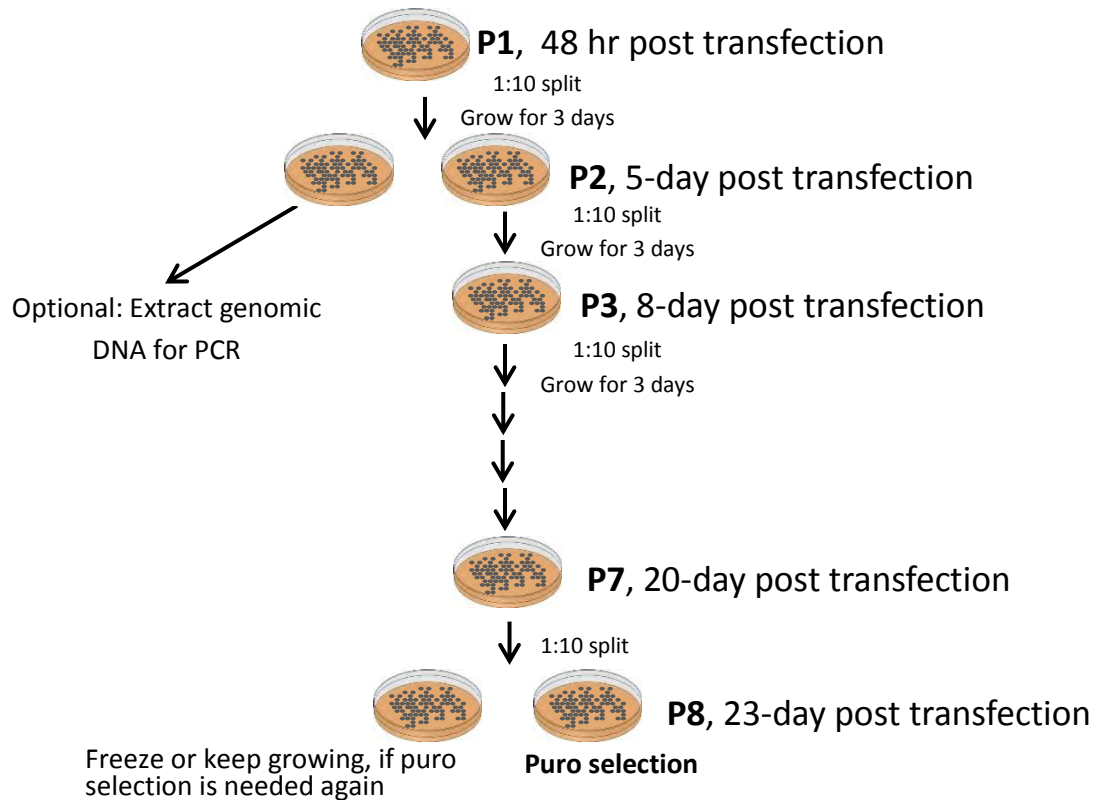


Figure 4. Diagram of genomic PCR Primer design.

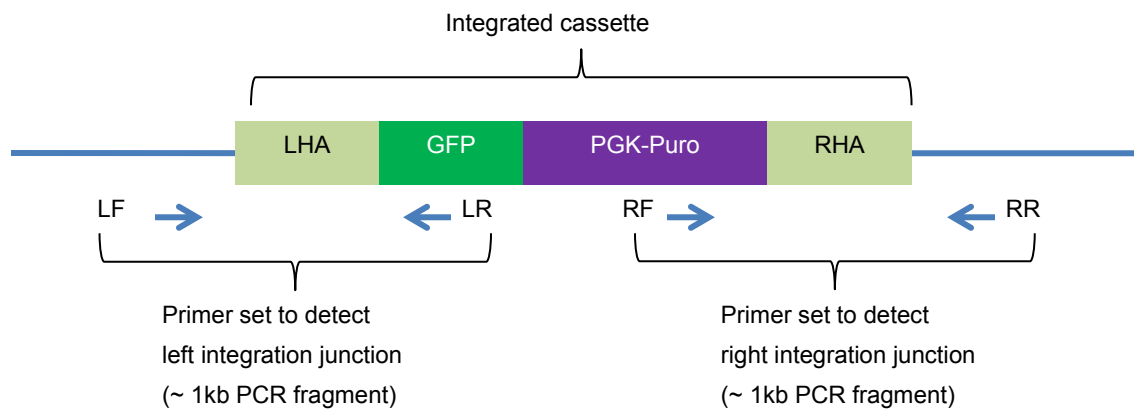

LF, LR: Forward and reverse PCR primer to amplify the left integration junction  
 RF, RR: Forward and reverse PCR primer to amplify the right integration junction

## **KN2.0, non-homology mediated CRISPR knockout kits**

### **Package contents**

- 2 vials of gRNA vectors, (SKU KNxxxxxxG1, KNxxxxxxG2), 3-5 µg DNA in TE buffer
- 1 vial of linear donor LoxP-EF1a-GFP-P2A-Puro-LoxP (SKU KNxxxxxD), 10 µg, lyophilized.
- Certificate of Analysis
- Application Guide: <https://www.origene.com/support/product-support/product-manuals>

*Note: The product is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, it is guaranteed to be stable for 12 months.*

### **Related Optional Reagents**

- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)

### **Related OriGene Products**

- Transfection reagent: <https://www.origene.com/products/others/transfection-reagents>
- Reagents and supplies for immunoblots: user preferred. OriGene has a selection of antibodies and detection reagents that are available at <https://www.origene.com/products/antibodies>
- DNA purification reagents <https://www.origene.com/products/others/nucleic-acid-purification-kits>
- qPCR reagents <https://www.origene.com/products/gene-expression/qpcr>
- CRISPR/Cas9 products <https://www.origene.com/products/gene-expression/crispr-cas9>


### **Notice to purchaser**

This product is for research use only. Use in and/or for diagnostics and therapeutics is strictly prohibited. By opening and using the product, the purchaser agrees to the following: The plasmids may not be distributed, resold, modified for resale or used to manufacture commercial products without prior written approval from OriGene Technologies, Inc. If you do not agree to the above conditions, please return the UNOPENED product to OriGene Technologies, Inc. within ten (10) days of receipt for a full refund.

### **Product Description**

KN2.0 gene knockout kit is non-homology mediated CRISPR kit. Target specific gRNA will cut the genome, then the linear donor DNA containing a selection cassette will be integrated at the gRNA cutting site at forward or reverse direction. The knockout efficiency is higher than HDR-mediated gene knockout. The majority gene knockout is bioallelic, one allele has donor integration, the other allele has indels (insertion and deletion). Indels might affect protein coding or cause premature stop.

Figure 5. Diagram of KN2.0 non-homology-mediated CRISPR knockout kit

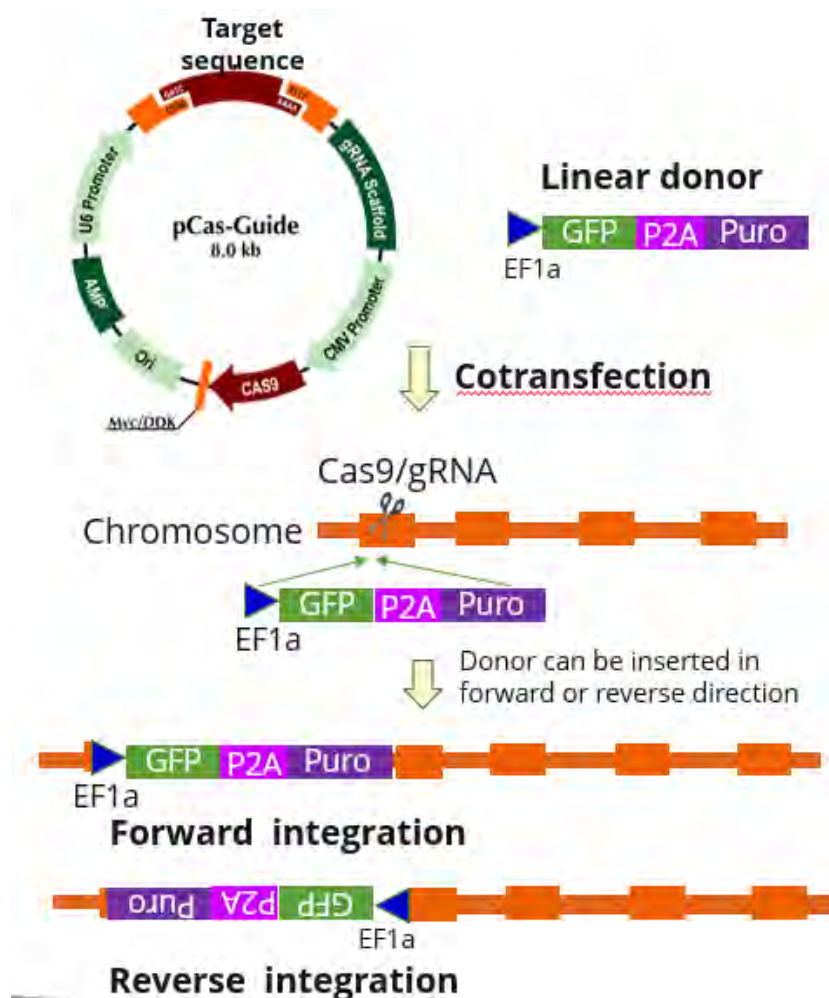

### Experimental Protocol

Each kit contains two gRNA vectors, and one donor DNA. To ensure high efficiency of cleavage, two gRNA constructs are provided.

A sample protocol listed below is for experiments performed in 6-well plates and using [TurboFectin](#) (cat# TF81001) as transfection reagent. If your experiments require a different size of culture plates, just scale up or down the reagents accordingly based on the relative surface area of your plate. Different type of cells may need a different transfection reagent; please follow the manufacturer's corresponding protocol. OriGene's Virus-like [Viromers](#) which are best for difficult-to-transfect cells.

1. Approximately 18-24 hours before transfection, plate  $\sim 3 \times 10^5$  adherent cells in 2 ml culture medium into each well of a 6-well plate or  $\sim 5 \times 10^5$  suspension cells per well to obtain 50-70% confluence on the following day. The amount of cells varies depending on the size of your cells.

2. Transfection in complete culture medium. Two separate transfections:



In a small sterile tube, combine the following reagents in the prescribed order. The order of reagent addition is important to achieve the optimal results.

- Dilute 1  $\mu\text{g}$  of one of the gRNA vectors in 250  $\mu\text{L}$  of Opti-MEM I (Life Technologies), vortex gently. Then add 1  $\mu\text{g}$  of the donor DNA into the same 250  $\mu\text{L}$  of Opti-MEM I. Vortex gently.
- Add 6  $\mu\text{L}$  of Turbofectin 8.0 to the diluted DNA (not the reverse order) and pipette gently to mix completely.
- Incubate the mixture 15 minutes at room temperature.

**Note:** 1. We recommend starting with the ratios of 3:1, Turbofectin:DNA. however, subsequent optimization may further increase the transfection efficiency.

**Note:** 2. For your specific cells, you may need to find the best transfection method and follow the manufacturer's protocol. Viromer transfection reagents work better for hard-to-transfect cells. In general, transfection efficiency of linear DNA is lower than circular DNA.

- Add the mixture above drop-wise to the cells plated in step 1 (no need to change the media). Gently rock the plate back-and-forth and side-to-side to distribute the complex evenly.
- Incubate the cells in a 5%  $\text{CO}_2$  incubator.

3. 48 hrs post transfection, split cells 1:10, grow additional 3 days; then split the cells again 1:10. Split cells 2-4 times in total. Since puromycin resistant gene in the donor DNA is under EF1a-P2A, the linear donor DNA before genomic integration will also provide puromycin resistance. The reason to grow cells for around 2 weeks before puromycin selection is to dilute out cells containing non-integrated donor DNA. The shorter cell passaging time before puromycin selection, the higher false positive rate. If you passage cells for 2 weeks after transfection, the correct gene knockout rate after puromycin selection is around 50%; if you only passage cells for 1 week, the gene knockout rate is around 30%.

#### Time lines of genome editing

- ✓ CRISPR targeted gene knockout / knockin--- 1 week post transfection
- ✓ Non-integrated donor DNA dilution with cell passaging--- 2 weeks post transfection

**Note 1.** Since stable cell selection takes time, you can try to analyze the cells at P2 (passage 2) to detect genomic integration using genomic PCR (Fig 6). Primer pair needs to be designed to

amplify the donor-inserted allele due to the small percentage of edited cells. Therefore, one primer should be outside of the donor cassette and one primer is in the donor cassette. There could be some difficulties in the genomic PCR at this step before selection due to the small percentage of edited cells and difficulties of genomic PCR. qPCR measuring the targeted mRNA level would not work due to the small percentage of edited cells.

4. Apply puromycin selection or GFP sorting. Split P3 or P5 cells 1:10, then grow cells directly in the puromycin containing complete media in 10 cm dishes (apply puro selection at P6 might work, but you might get more random donor integrated clones than P8). The dose of puromycin needs to be determined with a kill curve to find out the lowest dose that kills the non-transfected cells completely 4-7 days post selection; the range of puromycin is 1µg/ml to 10µg/ml). Change the media every 2-3 days.

**Note:** We recommend you still keep growing or freeze some of the transfected cells without selection; just in case, you need to perform the selection again.

5. The puromycin resistant cells are ready to be analyzed for genome editing. WB can be used to measure gene knockdown if there is a good protein-specific antibody available (better after isolating individual cell colonies). You will also need to do genomic PCR to verify the integration of the functional cassette (figure 6). With primer pair of 5F and 3R, both alleles of donor inserted and non-edited/indel will be amplified. Please see the data in figure 7. The smaller PCR fragment could be un-edited allele or allele containing indels that could cause protein reading frame change or premature stop. You can directly sequence the amplified genomic fragment using the PCR primers to verify the sequence. It is better after isolating single cell colonies, so the sequencing data is clear.


Figure 6. Diagram of PCR primers for genomic PCR to verify donor insertion.

### Forward integration



To detect 5' junction, use primer pairs 5F & 5R,  
3' junction, use primer pairs 3F & 3R

### Reverse integration

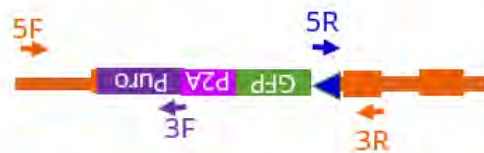

To detect 5' junction, use primer pairs 5F & 3F,  
3' junction, use primer pairs 5R & 3R

6. Isolate individual cell colonies.

There are two main methods to isolate single cell colonies, limiting dilution and cloning rings / cylinder.

#### 1) Limiting dilution.

This method is better to be used after puromycin selection. Dilute cells to seed about 1-2 cells/well in 96-well plate, after 1-2 weeks, observe under the microscope and select the wells only containing one cell colony, then further expand them to 6-well plate when they are confluent in the 96-well plate and so




on.

2) Cloning rings / cylinder

this method can be used in the same time with puromycin selection. Seed cells at lower density, such as 5% confluency in a larger cell culture dish, such as 10cm dish, apply puromycin selection when seeding.

**Note.** After puromycin selection, we observed some colonies are green, some colonies are not green. The efficiency of donor insertion is similar for the green and dark colonies. The reasons that the GFP expression is low in some cell colonies are not clear. It is ok to use GFP to sort out edited cells, so you will save one week comparing to puromycin selection. When you use puromycin selection, we recommend you analyze both green and dark cell colonies.

Figure 7. Genomic PCR verification data using primer pair 5F and 3R to amplify both donor-inserted and non-edited/indel alleles.

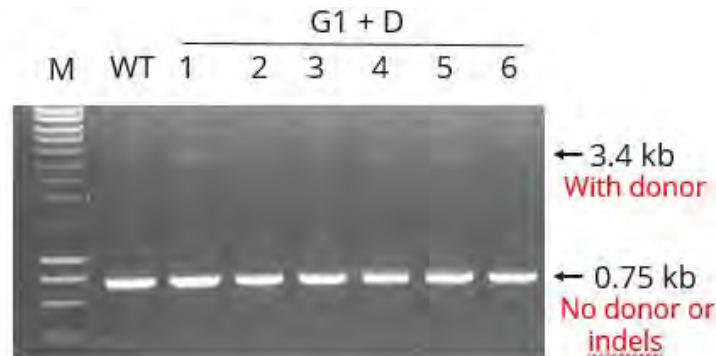

Single HEK293T cell colonies were isolated after puromycin selection. Genomic DNA was extracted and PCR was performed using primer pair 5F and 3R. WT: untransfected cells. 1, 2, 3, 4, 5, 6: single colonies of G1 and donor transfected.

**EF1a-GFP-P2A-Puro selection cassette sequence:**

ATAACTTCGTATAATGTATGCTATACGAAGTTATCGTGAGGCTCCGGTGCCCGTCAGTGGGCAGAGCGCACATCGC  
 CCACAGTCCCCGAGAAGTTGGGGGGAGGGGTGCGCAATTGAACCGGTGCCTAGAGAAGGTGGCGCGGGGTAAA  
 CTGGGAAAGTGATGTCGTGTACTGGCTCCGCCTTTTCCCGAGGGTGGGGGAGAACCGTATATAAGTGCAAGT  
 CGCCGTGAACGTTCTTTTTCGCAACGGGTTTGCCGCCAGAACACAGGTAAGTGCCGTGTGTGGTTCCCGCGGGCCT  
 GGCCTCTTACGGGTTATGGCCCTTGCCTGCTTGAATTACTTCCACCTGGCTGCAGTACGTGATTCTTGATCCCGA  
 GCTTCGGGTTGGAAGTGGGTGGGAGAGTTCGAGGCCTTGCCTTAAGGAGCCCCTTCGCCTCGTGCTTGAGTTGA  
 GGCCTGGCCTGGGCGCTGGGGCCCGCGCTGCGAATCTGGTGGCACCTTCGCGCCTGTCTCGCTGCTTTCGATAA  
 GTCTCTAGCCATTTAAAATTTTTGATGACCTGCTGCGACGCTTTTTTCTGGCAAGATAGTCTTGAAATGCGGGCC  
 AAGATCTGCACACTGGTATTTTCGGTTTTTGGGGCCGCGGGCGGCGACGGGGCCCGTGCCTCCAGCGCACATGTT  
 CGGCGAGGCGGGGCTGCGAGCGCGCCACCGAGAATCGGACGGGGTAGTCTCAAGCTGGCCGGCCTGCTCT  
 GGTGCTGGCCTCGCGCCCGCTGTATCGCCCCGCCCTGGGCGGCAAGGCTGGCCCCGTCGGCACCAAGTTGCGTG  
 AGCGGAAAGATGGCCGCTTCCCGGCCCTGCTGCAGGGAGCTCAAATGGAGGACGCGGCGCTCGGGAGAGCGG  
 GCGGGTGAGTACCCACACAAAGGAAAAGGGCCTTCCGTCCTCAGCCGTCGCTTCATGTGACTCCACGGAGTAC

CGGGCGCCGTCCAGGCACCTCGATTAGTTCTCGAGCTTTTGGAGTACGTCTTTAGGTTGGGGGAGGGGTTT  
TATGCGATGGAGTTTCCCCACACTGAGTGGGTGGAGACTGAAGTTAGGCCAGCTTGGCACTTGATGTAATTCTCT  
TGGAATTTGCCCTTTTTGAGTTTGGATCTTGGTTCACTCAAGCCTCAGACAGTGGTTCAAAGTTTTTTCTTCCAT  
TTCAGGTGTCGTGAATGGAGAGCGACGAGAGCGGCCTGCCCGCCATGGAGATCGAGTGCCGCATCACCGGCACC  
CTGAACGGCGTGGAGTTCGAGCTGGTGGGCGGCGGAGAGGGCACCCCGAGCAGGGCCGCATGACCAACAAGA  
TGAAGAGCACCAAAGGCGCCCTGACCTTCAGCCCCTACCTGCTGAGCCACGTGATGGGCTACGGCTTCTACCACTT  
CGGCACCTACCCAGCGGCTACGAGAACCCTTCTGCACGCCATCAACAACGGCGGCTACACCAACACCCGCATC  
GAGAAGTACGAGGACGGCGGCGTGTGCACGTGAGCTTACGCTACCGCTACGAGGCCGCGCCGCTGATCGGCGA  
CTTCAAGGTGATGGGCACCGGCTTCCCCGAGGACAGCGTGATCTTACCGACAAGATCATCCGCAGCAACGCCAC  
CGTGGAGCACCTGCACCCATGGGCGATAACGATCTGGATGGCAGCTTACCCGCACCTTACGCCTGCGCGACGG  
CGGCTACTACAGCTCCGTGGTGGACAGCCACATGCACTTCAAGAGCGCCATCCACCCAGCATCTGCAGAACGG  
GGCCCCATGTTCCGCTTCCGCCGCGTGGAGGAGGATCACAGCAACACCGAGCTGGGCATCGTGGAGTACCAGC  
ACGCCTTCAAGACCCCGGATGCAGATGCCGGTGAAGAAAGAGGAAGCGGAGCTACTAACTTCAGCCTGCTGAAG  
CAGGCTGGAGACGTGGAGGAGAACCCTGGACCTATGACCGAGTACAAGCCCACGGTGCGCCTCGCCACCCGCGA  
CGACGTCCCAGGGCCGTACGCACCCTCGCCGCGGTTCCGCCACTACCCGCCACGCGCCACACCGTTCGATCCG  
GACCGCCACATCGAGCGGGTACCGAGCTGCAAGAACTTCTCTCACGCGCGTCCGGCTCGACATCGGCAAGGTG  
TGGGTGCGGACGACGGCGCCGCGGTGGCGGTCTGGACCACGCCGAGAGCGTGAAGCGGGGGCGGTGTTCCG  
CCGAGATCGGCCCGCATGGCCGAGTTGAGCGGTTCCCGGCTGGCCGCGCAGCAACAGATGGAAGGCCTCCTG  
GCGCCGACCCGGCCAAAGGAGCCCGGTGGTCTTGGCCACCGTCCGGCTCTCGCCGACCACAGGGCAAGGG  
TCTGGGACGCGCGTCTGCTCCCCGAGTGGAGGCGGCCGAGCGCGCCGGGTGCCCGCTTCTGGAGACCT  
CCGCGCCCCGCAACCTCCCCTTCTACGAGCGGCTCGGCTTACCGTACCGCCGACGTCGAGGTGCCGAAGGACC  
GCGCACCTGGTGCATGACCCGCAAGCCCGGTGCCTGAACTGTTTATTGCAGCTTATAATGGTTACAAATAAAGC  
AATAGCATCACAAATTTACAAATAAAGCATTTTTTTCACTGCATTCTAGTTGTGGTTTGTCCAAACTCATCAATGTA  
TCTTAATAACTTCGTATAATGTATGCTATACGAAGTTAT

LoxP in cyan, EF1A promoter in blue, tGFP (699 bp) in green, Puromycin resistance gene (600 bp) purple. P2A sequence (66 bp) is in red.

## II. Precut All-in-one CRISPR vector cloning kits (pCas-Guide (SKU: GE100001), pLenti-Cas-Guide (SKU:GE100009))

### Package contents

The following components are included:

- One (1) vial of precut CRISPR vector DNA, lyophilized ready for ligation (10 RXNs). Reconstitute in 10  $\mu$ L dH<sub>2</sub>O, final concentration 10 ng/  $\mu$ L.
- One (1) vial of forward sequencing primer, CF3 (SKU GE100008) to sequence the targeting sequence cloned into the CRISPR vector; 100 picomoles dried onto the bottom of screw cap tubes. Reconstitute the primer in 10  $\mu$ L dH<sub>2</sub>O to make a 10  $\mu$ M solution.
- One (1) vial of annealing buffer (SKU GE100007) for oligo annealing (10X), 200  $\mu$ L
- Certificate of Analysis
- Application Guide

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit <https://www.origene.com/products/others/nucleic-acid-purification-kits> ). The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

### **Related Optional Reagents**

Nuclease free water  
 T4 DNA ligase and buffer  
 Competent *E. coli* cells  
 LB agar plates with ampicillin, 100 µg/mL  
 LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)  
 DNA purification reagents

### **Related OriGene Products**

- TrueClone™ FL cDNA clones <https://www.origene.com/products/cdna-clones/trueclone>
- HuSH™ shRNA Plasmids <https://www.origene.com/products/rnai/shrna-plasmids>
- Validated Antibodies <https://www.origene.com/products/antibodies>
- Purified Proteins <https://www.origene.com/products/proteins>
- Over-expression lysates <https://www.origene.com/products/proteins/over-expression-lysates>
- Transfection Reagents <https://www.origene.com/products/others/transfection-reagents>
- Anti-tag Antibodies <https://www.origene.com/products/antibodies/tag-antibodies>


### **Notice to purchaser**

This product is for research use only. Use in and/or for diagnostics and therapeutics is strictly prohibited. By opening and using the product, the purchaser agrees to the following: The plasmids may not be distributed, resold, modified for resale or used to manufacture commercial products without prior written approval from OriGene Technologies, Inc. If you do not agree to the above conditions, please return the UNOPENED product to OriGene Technologies, Inc. within ten (10) days of receipt for a full refund by contacting customer service at [custsupport@origene.com](mailto:custsupport@origene.com).

### **Production and Quality Assurance**

The precut pCas-Guide or pLenti-Guide vectors have been tested to successfully religate to annealed oligo DNA fragments. When OriGene experimental protocol is followed (details on page 13-15), 1 µL of the ligation reaction generated with the precut vector can produce 100 colonies when transformed into 10<sup>6</sup> cfu/µg competent cells. The self-ligation background (vector religation to itself without an insert) is less than 5% of transformants. The amount of digested DNA provided in the kit is sufficient for ten ligation reactions.

Figure 5. The vector map of pre-cut pCas-Guide and pLenti-Cas-Guide

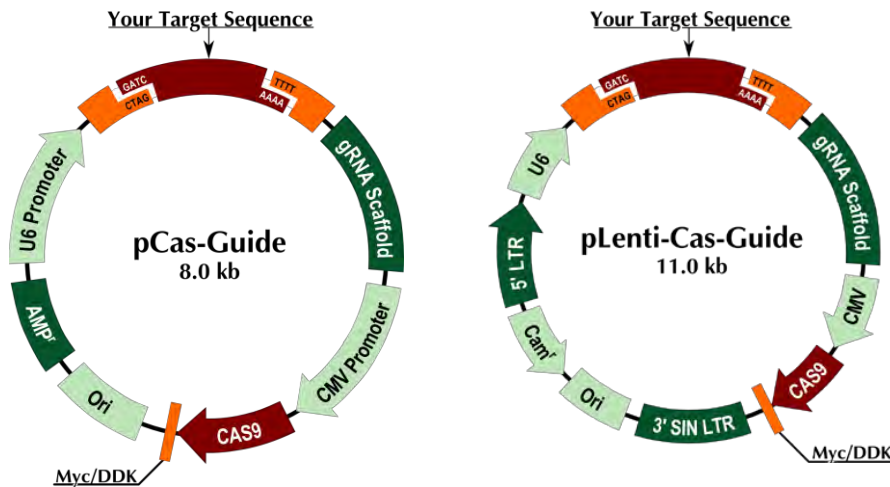

### Product Description

Both pCas-Guide and pLenti-Cas-Guide vectors are designed for cloning a guide RNA insert for genome editing purpose. The vectors also have a CMV-driven codon-optimized Cas9. When co-transfected with a proper donor DNA, targeted genome editing can be achieved. The vectors are supplied as pre-cut vectors, ready for insert ligation. This system has been successfully validated in multiple cases of genome editing.

### Experimental Protocols

#### I. Design target sequence

OriGene has developed a proprietary Cas9 guide RNA designing tool which is free to use, <http://www.blueheronbio.com/>. Follow the instructions below to design your guide RNA:

1. Select your desired Cas9 cutting site from your genomic region of interest.
2. Copy around 100 bp of genomic sequence flanking the cutting site (-50 to +50). Paste the sequence to the sequence box and click the Search button.
3. The program will return all possible targeting sequences with location and GC content obtained from searching both the plus and minus strands. If there is no target returned, expand your genomic region of interest (-100 to +100) and search again until there is a positive return.
4. Select a few target sequences to Blast against the genomic DNA database to check sequence specificity.
5. Select 2 to 3 target sequences to clone into pCas-Guide vector.

#### II. Addition of extra bases to the ends of the target sequence

To facilitate cloning of the 20-bp target sequence, extra bases need to be added to the ends.

1. Select a desired 20-bp sequence as a target. The following is an example sequence:

Forward sequence: 5' ATGGGAGGTGGTATGGGAGG 3'  
 Reverse complement sequence: 5' CCTCCCATACCACCTCCCAT 3'

2. Add 'gatcg' to the 5' end of the forward sequence and 'g' to its 3' end. The final sense oligo in this example will be

5' gatcgATGGGAGGTGGTATGGGAGGg 3'

3. Add 'aaaac' to the 5' end of reverse complementary sequence and 'c' to its 3' end. The final reverse complementary sequence is

5' aaaacCCTCCCATACCACCTCCCATc 3'

The two oligos should anneal to form the following double strand:

5' gatcgXXXXXXXXXXXXXXXXXXXXXXXXXg 3'  
 3' cXXXXXXXXXXXXXXXXXXXXXXXXXcaaaa 5'

4. Order the two final oligos from a commercial oligo provider, such as IDT.

### III. Cloning the double-stranded oligos into the pCas-Guide vector

1. Anneal the oligos to form double-stranded duplexes

In a PCR tube, add the following:

2 µL Forward oligo (100 µM stock)

2 µL Reverse oligo (100 µM stock)

4 µL 10X annealing buffer

32 µL dH<sub>2</sub>O

Mix the solution and follow the steps to anneal the oligos in a PCR machine:

94°C for 4min

75°C for 5 min

65°C for 15 min

25°C for 20 min

After annealing, transfer the solution to a 1.5 mL tube and add 360 µL of dH<sub>2</sub>O.

The double-stranded oligo DNA is ready for ligation.

2. Ligation and transformation

A. Prepare the ligation according to the following protocol

Component	Volume
10x Ligation buffer	1 µL
Precut pCAS-Guide vector (10 ng/ µL)	1 µL
Annealed double-stranded oligos (diluted from step 1)	1 µL
Ligase (0.5 u/ µL, Weiss unit)	0.5 µL
dH <sub>2</sub> O	6.5 µL
<b>Total Volume</b>	<b>10 µL</b>

- B. Mix the solution and incubate the tube at 22 to 37°C or room temperature for two hours according to the manufacturer's recommendation.
- C. Add 1 µL of the ligation mixture to 10 µL of competent cells (efficiency rated > 10<sup>6</sup> cfu/µg DNA) on ice. Do the transformation according to the manufacturer's protocol. For chemically competent cells, follow steps D-E.
- D. Mix the tube gently and keep it on ice for 25 minutes.
- E. Heat shock the tube for 30 seconds at 42°C.
- F. Put the tube on ice for 2 minutes, then add 500 µL LB or SOC medium.
- G. Rock the tube gently at 37°C for 1 hour.
- H. Spread 50 µL of the *E. Coli* cells on an LB agar plate containing the corresponding Antibiotics (most vectors are ampicillin resistant except lenti CRISPR vectors, which is chloramphenicol).
- I. Centrifuge the remaining *E. Coli* cells at 5K rpm for 5 minutes. Discard the majority of the supernatant (around 50 µL supernatant left) and resuspend the cell pellet in the remaining liquid. Spread all the *E. Coli* cells on a separate LB-agar plate.
- J. Incubate the two plates at 37°C for 16 hours to allow colony formation.

### 3. Screening colonies

In a typical subcloning ligation, at least 95% of the colonies should contain the desired insert. Pick 6 to 10 colonies into 5 mL LB-ampicillin culture each, and culture overnight. Perform DNA purification using a mini-prep kit from OriGene, <https://www.origene.com/products/others/nucleic-acid-purification-kits/>. Sequence the purified DNA and analyze the sequencing data to identify a correct clone for proper insert identification and orientation.

#### **Lenti-based protocols:**

NOTE: Performing Lentiviral experiments REQUIRES special laboratory conditions and/or permissions (BL2). Follow the guidelines and regulations of your institution. Perform the experiments with due caution to avoid exposure to infectious materials.

- A. Production of pseudovirus (10 cm plate format, the production size can be scaled up or down accordingly):
  - 1. Day 1, plate HEK293T cells in a 10 cm dish to approximately 40% confluency the day before transfection (antibiotic-free preferred). Cells should reach 65-70% confluency within 24 hours.
  - 2. Day 2, In sterile tube, dilute the following DNA in 1.5 mL Opti-MEM, and pipet gently to mix completely.

- a. 5 µg of lenti vector, Lenti CRISPR vector or Lenti-ORF or lenti-shRNA
- b. 6 µg of packaging plasmids, cat# [TR30037](#) (for regular lentivirus production) or [TR30036](#) (for integration-deficient lentivirus production).

Add 33 µL of TurboFectin transfection reagent to the diluted DNA (not the reversed order), pipet gently to mix completely.

- c. Incubate for 15 min at room temperature.
- d. Add the transfection mixture prepared above dropwise to the cells. Gently rock the plate back-and-forth and from side-to-side to distribute the complex evenly. Incubate at 37° in a CO2 incubator.

Note: With TurboFectin, no medium change is necessary, directly add the transfection mixture to cells in complete growth media.

3. Day 3, change the growth medium and continue to incubate the plate for 48 hours.
4. Day 5
  - a. After the 48 hour incubation, transfer the cell culture supernatant to a 15 mL centrifuge tube.
  - b. Centrifuge the tubes at 3K RPM for 10 mins and filter the supernatant through a syringe filter (0.45 micron) and collect the viral solution to a new sterile tube.
5. The viral particles are ready to be used. If virus titer is needed, you can use PCR or Elisa methods and follow the corresponding protocols. Lentiviral particles can be stored at 4 °C for 2 weeks or store at -80 °C for long-term storage.

#### B. Transduction of lentivirus to target cells

1. Day 1, plate target cells in three 10 cm plates at a density that will produce approximately 60% confluency in 24 hours. Note: other size formats can also be used depending on the nature of your experiment. Adjust the reagent amount accordingly.
2. Day 2, Remove the growth media from the plates prepared the day before. To plate 1, add 4.5 mL of fresh growth medium and 0.5 mL of Lentiviral particles; To plate 2, add 4.0 mL of growth medium and 1 mL of Lentiviral particles; To plate 3, add 2.5 mL of growth medium and 2.5 mL of Lentiviral particles (for a low titer viral preparation, the amount of virus added can be increased to 5 mL). Mix the solution by gentle swirling.
3. Add 5 µl polybrene (1,000x, 8 mg/mL) to each plate. Mix by gentle swirling.
4. Incubate the cells at 37 °C with 5% CO2 for 4 hours. Remove the transduction medium and add 10 mL of fresh growth medium. Incubate the cells for three more days.

The transduced cells are ready for downstream analyses such as RNA and protein detection.

### III. Circular All-in-one CRISPR vectors (SKU: GE100002, GE100018, GE100022, GE100010, GE100045)

#### Package contents

- One (1) vial of CRISPR vector DNA, 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- One (1) vial of CF3 sequencing primer (SKU GE100008) is provided for vectors (GE100002, GE100010) to sequence the targeting sequence cloned into CRISPR vector; 100 picomoles dried onto the bottom of screw cap tubes. Reconstitute the primer in 10 µL dH<sub>2</sub>O to make a 10 µM solution.
- Certificate of Analysis
- Application Guide available online.

*\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit*

*<https://www.origene.com/products/others/nucleic-acid-purification-kits> )*

*The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.*

#### Related Optional Reagents

- Oligo annealing buffer, SKU [GE100007](#)
- BamH I and BsmB I
- T4 DNA ligase and buffer
- Competent *E. coli* cells, <https://www.origene.com/products/others/competent-cells>
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

#### Related OriGene Products

- Cloning the homologous arms in predesigned donor rescue construct harboring functional cassettes (service from [www.blueheronbio.com](http://www.blueheronbio.com) ).
  - GFP-PGK-Loxp-Puro-Loxp
  - RFP-PGK-Loxp-BSD-Loxp
  - Luciferase-PGK-Loxp-Puro-Loxp
  - mBFP-PGK-Loxp-Neo-Loxp
- TrueClone™ FL cDNA clones <https://www.origene.com/products/cdna-clones/trueclone>
- HuSH™ shRNA Plasmids <https://www.origene.com/products/rnai/shrna-plasmids>
- Validated Antibodies <https://www.origene.com/products/antibodies>
- Purified Proteins <https://www.origene.com/products/proteins>
- Over-expression lysates <https://www.origene.com/products/proteins/over-expression-lysates>
- Transfection Reagents <https://www.origene.com/products/others/transfection-reagents>
- Anti-tag Antibodies <https://www.origene.com/products/antibodies/tag-antibodies>

#### Notice to purchaser

This product is for research use only. Use in and/or for diagnostics and therapeutics is strictly prohibited. By opening and using the product, the purchaser agrees to the following: The





plasmids may not be distributed, resold, modified for resale or used to manufacture commercial products without prior written approval from OriGene Technologies, Inc. If you do not agree to the above conditions, please return the UNOPENED product to OriGene Technologies, Inc. within ten (10) days of receipt for a full refund by contacting customer service at [custsupport@origene.com](mailto:custsupport@origene.com).

### Product Description

The CRISPR vectors are designed for cloning a guide RNA insert for genome editing purpose. Target sequence can be cloned into the vector via BamH I and BsmB I sites. The vectors also express a CMV-driven codon-optimized Cas9. When co-transfected with a proper donor DNA, your desired targeted genome editing can be achieved. The vectors retain the ampicillin resistance gene for the selection of *E. coli* transformants.

Figure 6. The vector maps of the circular CRISPR vectors.

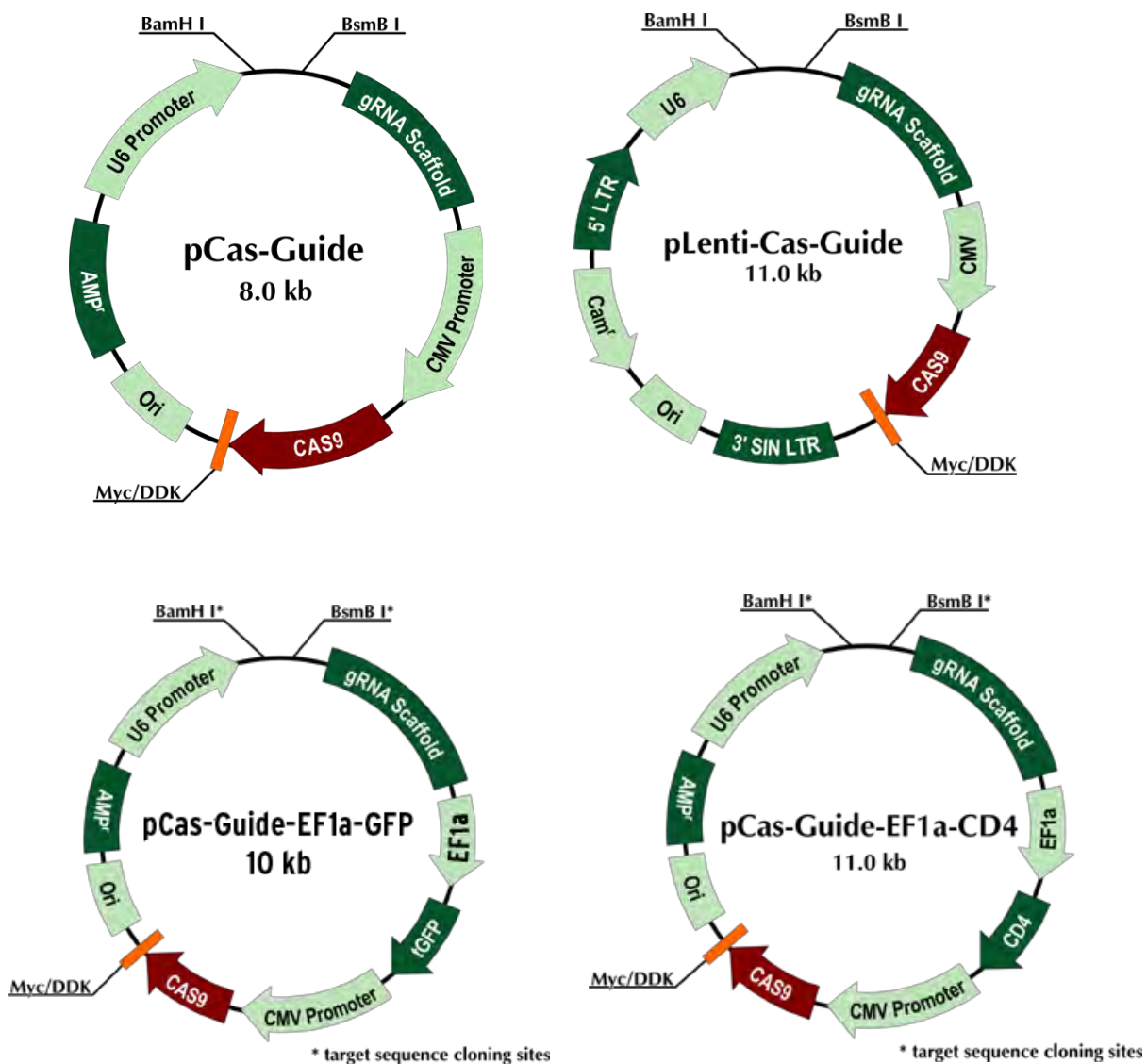

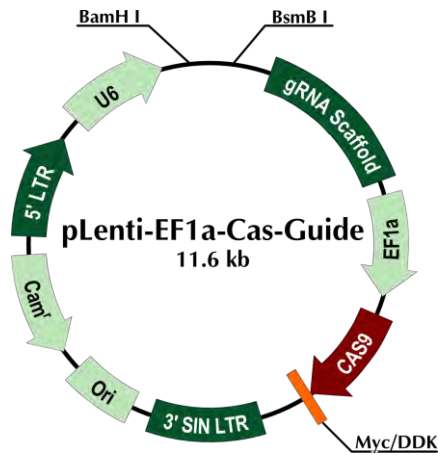

### Experimental protocol

1. Digest pCas-Guide plasmid with BamH I and BsmB I

Resuspend the 10 µg lyophilized DNA in 100 µL dH<sub>2</sub>O, and incubate for at least 30 min before use. Set up a digestion reaction as described below.

Component	Volume
10X restriction buffer	3 µL
BamH I	0.8 µL
BsmB I	0.8 µL
Nuclease free water	15.4 µL
Vector DNA	10 µL
Total volume	30 µL

Incubate the reaction at 37°C for 3 hrs, then add 1 µL antarctic phosphatase (units used according to the manufacturer's protocol), and continue the incubation at 37°C for another 30 min. Dephosphorylation of the digested vector is essential to eliminate self-ligation. Purify the desired vector fragment by running the digestion reaction on an agarose gel, and isolate the appropriate band using a gel purification column. Elute the digested plasmid vector in 40 µL of 10 mM Tris buffer.

2. Target sequence designing and cloning into the pre-cut vectors, please follow the detailed protocol from page 13-17 in this manual.

Lenti based protocol, please refer to page 16-17.

#### IV. gRNA only Vectors, pGuide (SKU GE100042), pGuide-EF1a-GFP (SKU GE100044), pLenti-Guide-Puro (SKU GE100032)

##### **Package contents**

- One (1) vial of circular pGuide or pLenti-Guide-Puro plasmid DNA, 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit

<https://www.origene.com/products/others/nucleic-acid-purification-kits> )

The cDNA clone is shipped at room temperature, but should be kept at -20°C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

##### **Related Optional Reagents**

- Nuclease free water
- Oligo annealing buffer, SKU [GE100007](#)
- BamH I, BsmB I
- T4 DNA ligase and buffer
- Competent *E. coli* cells, <https://www.origene.com/products/others/competent-cells>
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

##### **Related OriGene Products**


- Cloning the homologous arms in predesigned donor rescue construct harboring functional cassettes (service from <https://www.origene.com/products/gene-expression/crispr-cas9/grna-cloning-donor-construction> ).

GFP-PGK-Loxp-Puro-LoxP  
 RFP-PGK-Loxp-BSD-LoxP  
 Luciferase-PGK-Loxp-Puro-LoxP  
 mBFP-PGK-Loxp-Neo-LoxP

##### **Product Description**

gRNA only vectors can be used if your target cells already expressing Cas9. Two gRNA only cloning vectors are offered, pGuide, in regular mammalian expression vector, pLenti-Guide-Puro, in lenti vector backbone and contains puromycin selection marker. Target sequence can be cloned into the vector via BamH I and BsmB I sites. pGuide retains the ampicillin resistance gene for the selection of *E. coli* transformants; pLenti-Guide-Puro retains chloramphenicol resistance for *E. coli* selection.

Figure 7. Plasmid maps of pGuide and pLenti-Guide-Puro.

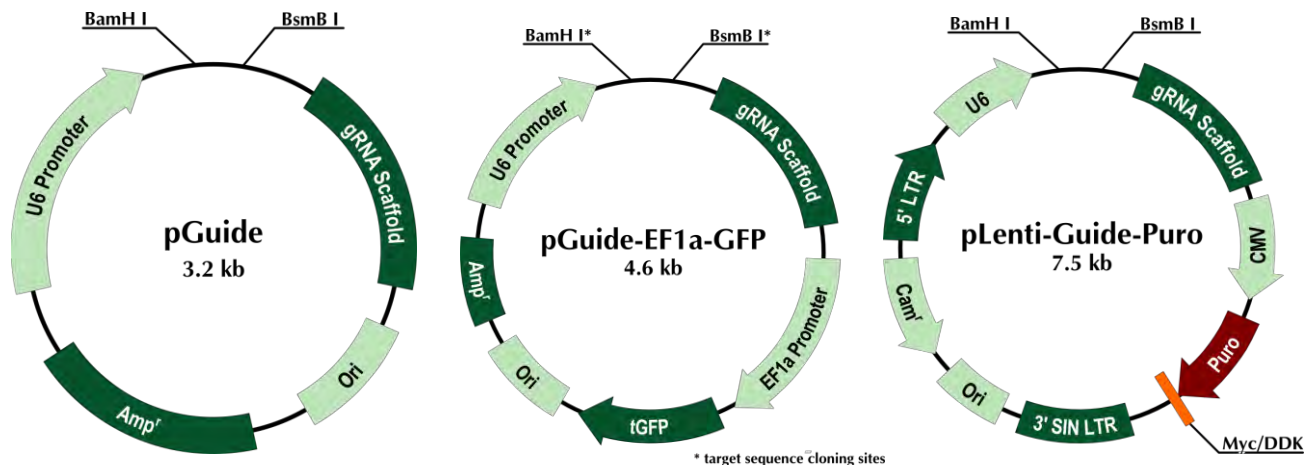

### Experimental protocol

The protocol of the gRNA only cloning is the same as the all-in-one CRISPR vector. Please see details from page 13-17 in this manual.

### V. Cas9 only vectors,

- pLenti-Cas9 (SKU GE100028),**
- pLenti-Cas9-IRES-Puro (SKU GE100029),**
- pLenti-EF1a-Cas9-IRES-Puro (SKU GE100030),**
- pLenti-Cas9-P2A-tGFP (SKU GE100031)**
- pAAVS1-Cas9-Puro-DNR (SKU GE100037)**
- pAAVS1-Cas9-BSD-DNR (SKU GE100039)**

### Package contents

- One (1) vial of circular plasmid vector DNA, 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit

<https://www.origene.com/products/others/nucleic-acid-purification-kits>

The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

### Related Optional Reagents

- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

### Related OriGene Products

- gRNA only vectors, pGuide (GE100042), pLenti-Guide-Puro (SKU GE100032)
- Cloning the homologous arms in predesigned donor rescue construct harboring functional cassettes (service from <https://www.origene.com/products/gene-expression/crispr-cas9/grna-cloning-donor-construction> ).
  - GFP-PGK-Loxp-Puro-LoxP
  - RFP-PGK-Loxp-BSD-LoxP
  - Luciferase-PGK-Loxp-Puro-LoxP
  - mBFP-PGK-Loxp-Neo-Loxp


### Product Description

#### Lenti Cas9 vectors

The Lenti Cas9 only vectors express Cas9 after being transfected (using transfection reagent) or transduced (packaged into lentiviral particles) into cells, no gRNA expression cassette in the vector. The Lenti Cas9 vectors can be transfected or transduced with gRNA construct into cells or animals at the same time or develop Cas9 stable cells first, then introduce sgRNA into cells. The vectors retain the chloramphenicol resistance gene for the selection of *E. coli* transformants.

Vector	pLenti-Cas9	pLenti-Cas9-IRES-Puro	pLenti-EF1a-Cas9-IRES-Puro	pLenti-Cas9-P2A-tGFP
Cas9 Promoter	CMV	CMV	EF1a	CMV
Mammalian selection	none	puro	puro	tGFP

Figure 8. Plasmid maps of Lenti Cas9 vectors

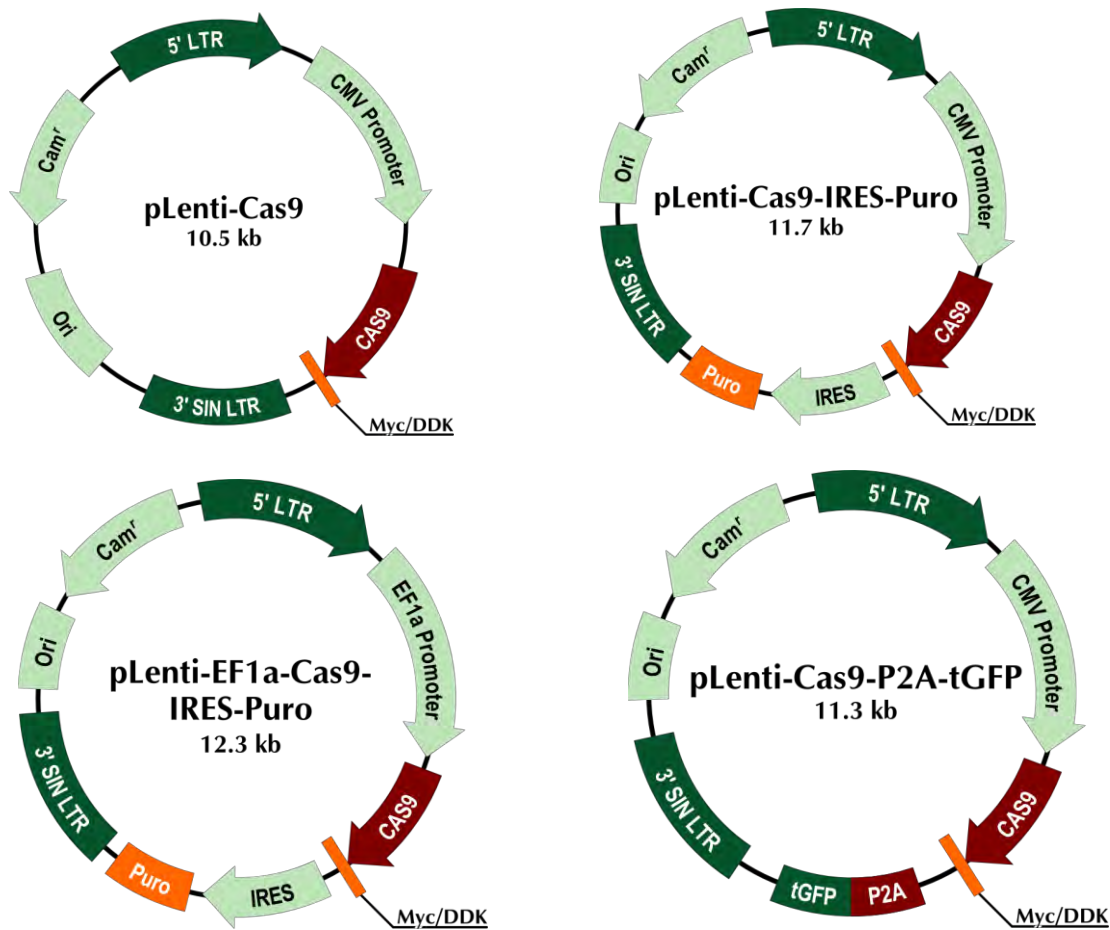


### pAAVS1-Cas9-Puro-DNR, pAAVS1-Cas9-BSD-DNR vectors

These two vectors can be used as regular mammalian Cas9 expression vectors when transfected into cells alone as Cas9 is driven by CMV promoter.

There are also AAVS1 homologous arm sequences in the vector, the vectors can be used as donor vectors together with pCas-guide-AAVS1 (gRNA targeting AAVS1, SKU GE100023) to insert Cas9 at AAVS1 locus in human cells. AAVS1 locus has been proven to be a safe harbor location for exogenous insertion to achieve stable and robust gene expression.

For details, please see page 44-45 in this manual.

Fig. 9. Diagram of how Cas9 is inserted at AAVS1 locus

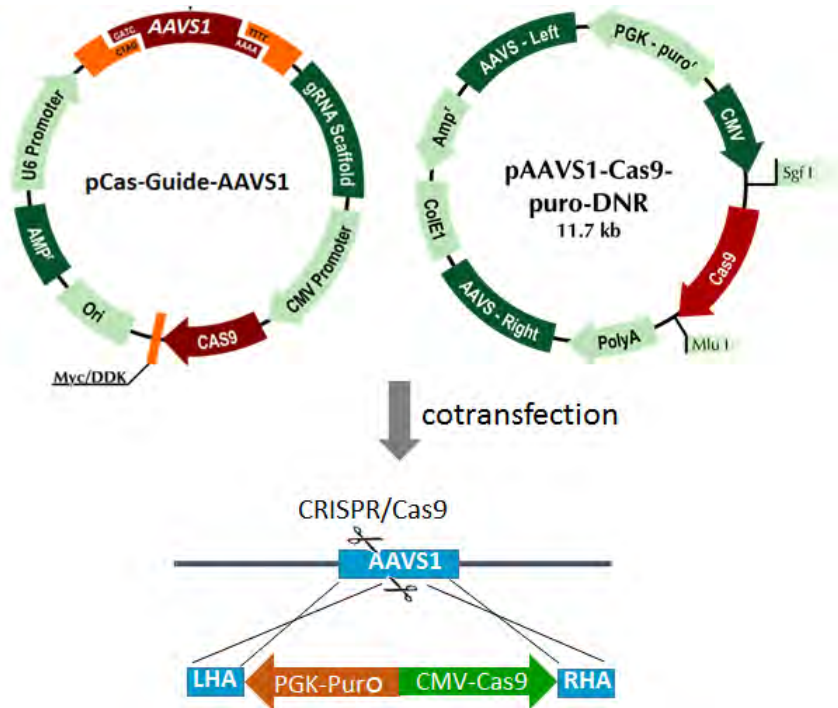

### Experimental Protocol for Lenti Cas9 Vectors

1. Transient transfection with gRNA construct. Follow the manual manufacture's protocol of the transfection reagents.
2. Transduction with gRNA particles. Making lentiviral particles and transducing cells, follow protocol on page 16-17 in this manual.

### VI. T7 driven CRISPR/Cas system-pT7-Guide-IVT (SKU: GE100025) and pT7-Cas9 (SKU: GE100014)

#### Package contents

- One (1) vial of circular pT7-Guide-IVT (SKU: GE100025) plasmid DNA or pT7-Cas9 (SKU: GE100014), 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit

<https://www.origene.com/products/others/nucleic-acid-purification-kits> )


The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

#### Related Optional Reagents

- Nuclease free water





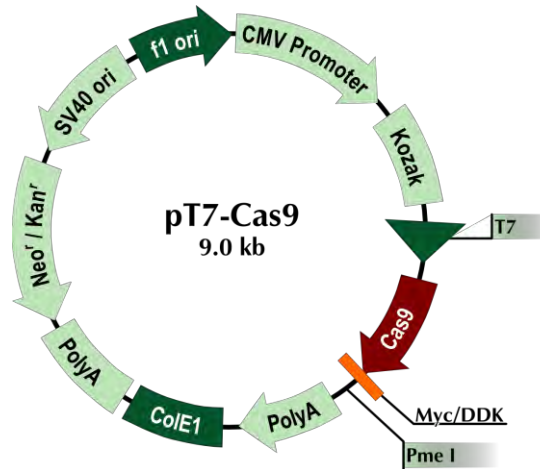

## Experimental protocol

### I. Design genomic target sequence and cloning into pT7-Guide-IVT vector

1. Digest pT7-Guide-IVT vector with BsmB I restriction enzyme

Resuspend the 10 µg lyophilized DNA in 100 µL dH<sub>2</sub>O, and incubate for at least 30 min before use. Set up a digestion reaction as described below.

Component	Volume
10X restriction buffer	3 µL
BsmB I*	0.8 µL
Nuclease free water	16.2 µL
Vector DNA	10 µL
Total volume	30 µL

Incubate the reaction at 37°C for 3 hrs (the isoschizomer ESP3 I from Thermo Scientific can be used). Purify the desired vector fragment by running the digestion reaction on an agarose gel, and isolate the appropriate band using a gel purification column. Elute the digested plasmid vector in 40 µL of 10 mM Tris buffer.

*Note: Please do not dephosphorylate the vector as oligos are not phosphorylated.*

2. Target sequence designing and cloning.

OriGene has developed a proprietary Cas9 guide RNA designing tool which is free to use, <http://www.blueheronbio.com/>. Design a target sequence of 20 bp,

To facilitate cloning of the 20-bp target sequence, extra bases need to be added to the ends.

- 1). Select a desired 20-bp sequence as a target. The following is an example sequence:

Forward sequence: 5' ATGGGAGGTGGTATGGGAGG 3'  
 Reverse complement sequence: 5' CCTCCCATACCCACCTCCCAT 3'

Add 'atagG' to the 5' end of the forward sequence and 'G' to its 3' end.  
 The final sense oligo in this example will be  
 5' atagG ATGGGAGGTGGTATGGGAGGg 3'

2). Add 'aaaac' to the 5' end of reverse complementary sequence and 'C' to its 3' end.

The final reverse example complementary oligo:  
 5' aaaacCCTCCCATACCCACCTCCCATc 3'

The two oligos should anneal to form the following double strand:

5' atagGxxxxxxxxxxxxxxxxxxxxxxxxxxG 3'  
 3' CxxxxxxxxxxxxxxxxxxxxxxxxxxCaaaa 5'

3). Order the two final oligos from a commercial oligo provider, such as IDT. The oligos are ready to be cloned into the BsmB I digested pT7-Guide-IVT vector. Following the oligo cloning procedure on page 15-16 of this manual.

3. Sequencing the cloned target sequence in pT7-Guide-IVT can be done by the common M13 forward primer: 5' CGCCAGGGTTTTCCAGTCACGAC 3'

## II. Producing gRNA and Cas9 mRNA using T7 in vitro transcription kits

To make gRNA using pT7-Guide-IVT after the genomic target sequence cloned, we recommend using MeGAscript T7 kit (Life Technologies) and follow the manufacturer's protocol. pT7-Guide-IVT vector can be linearized using EcoR I. EcoR I will cut T7 and gRNA out; but you don't need to purify the fragment. You only need to clean it using a PCR purification column. Then follow the MeGAscript T7 kit protocol to produce gRNA. To make Cas9 mRNA using pT7-Cas9, we recommend using mMESSAGE Mmachine T7 ULTRA kit (Life Technologies) and follow the manufacturer's protocol. pT7-Cas9 vector can be linearized using Pme I which is at the 3' end of Cas9 sequence. You can then clean up the Pme I digested reaction using a PCR purification column; then follow the mMESSAGE Mmachine T7 ULTRA kit protocol to produce capped and polyadenylated Cas9 mRNA.

## VII. Cre expression vector for Cre-Lox recombination, pCMV6-Entry Cre (SKU: GE100017)

### Package contents

- One (1) vial of circular pCMV6-Entry-Cre (SKU: GE100017) plasmid DNA 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis
- Application Guide available on line

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit <https://www.origene.com/products/others/nucleic-acid-purification-kits> )

The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

#### Related Optional Reagents

- Nuclease free water
- Competent *E. coli* cells, <https://www.origene.com/products/others/competent-cells>
- LB agar plates with Kanamycin, 25 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

#### Related OriGene Products


- Genome wide gene knockout/knockin kit via CRISPR, <https://www.origene.com/products/gene-expression/crispr-cas9/knockout-kits>
- Cloning the homologous arms in predesigned donor rescue construct harboring functional cassettes (service from <https://www.origene.com/products/gene-expression/crispr-cas9/grna-cloning-donor-construction> ).

GFP-PGK-Loxp-Puro-Loxp  
 RFP-PGK-Loxp-BSD-Loxp  
 Luciferase-PGK-Loxp-Puro-Loxp  
 BFP-PGK-Loxp-Neo-Loxp

#### Product Description

In our predesigned donor vector cassette, the PGK-puro cassette is flanked by two LoxP sites, which can be floxed out by Cre recombinase. pCMV6-Entry-Cre is a Cre recombinase mammalian expression vector in which Cre is under CMV promoter. Expressed Cre contains a C-terminal Myc-DDK tag (DDK is the same as Flag tag).

Figure 11. Plasmid map of pCMV6-Entry-Cre

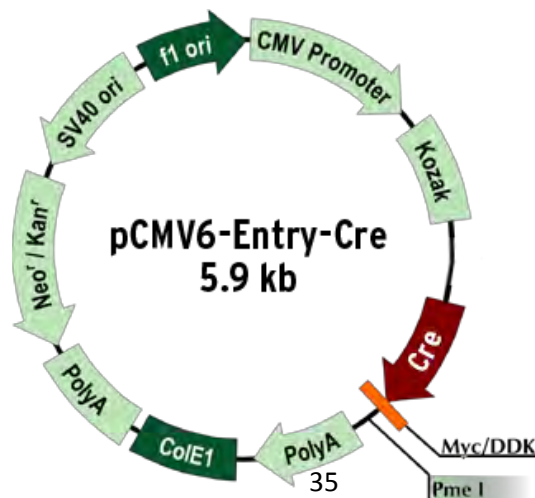

## VIII. Cas9 D10A nickase vectors, pCas-Guide-Nickase (SKU GE100019), pT7-Cas9-Nickase (SKU GE100020)

### Package contents

- One (1) vial of circular nickase plasmid DNA, pCas-Guide-Nickase (SKU: GE100019) or pT7-Cas9-Nickase, 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis
- Application Guide

### Related Optional Reagents

- Nuclease free water
- Oligo annealing buffer, SKU [GE100007](#)
- BamH I, BsmB I
- T4 DNA ligase and buffer
- [Competent E. coli cells](#)
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

### Related OriGene Products

- [Cloning the homologous arms](#) in predesigned donor rescue construct harboring functional cassettes.
  - GFP-PGK-Loxp-Puro-LoxP
  - RFP-PGK-Loxp-BSD-LoxP
  - Luciferase-PGK-Loxp-Puro-LoxP
  - mBFP-PGK-Loxp-Neo-LoxP


### Product description

WT Cas9 has two active nuclease domains and it can produce double-stranded DNA breaks. D10A mutation in Cas9 disables one nuclease domain; therefore Cas9D10A can only nick the targeted genome. pCas-Guide-Nickase (SKU GE100019) is in the same vector backbone as pCas-Guide (SKU GE100001 and GE100002) which is all-in-one vector, target sequence can be cloned and the vector express Cas9. pT7-Cas9D10A is in the same vector backbone as pT7-Cas9 (GE100014) which is used for *in vitro* production of Cas9 mRNA. Since Cas9D10A only nicks the genomic DNA, it needs two different gRNAs (one on sense strand and one on anti-sense strand) to cause double-stranded break. The off-target problem is significantly decreased as Cas9D10A needs two gRNAs.

### Experimental protocol

The experimental protocol of GE100019 will be similar to pCas-Guide, page 12 on this manual. The experimental protocol of GE100019 will be similar to pT7-Cas9, page 24 on this manual. Since Cas9D10A needs two different gRNAs, you will need to validate the cleavage efficiency of each gRNA using WT Cas9 before using them together with Cas9D10A vectors.

Figure 12. Plasmid maps of pCas-Guide-Nickase and pT7-Cas9-Nickase

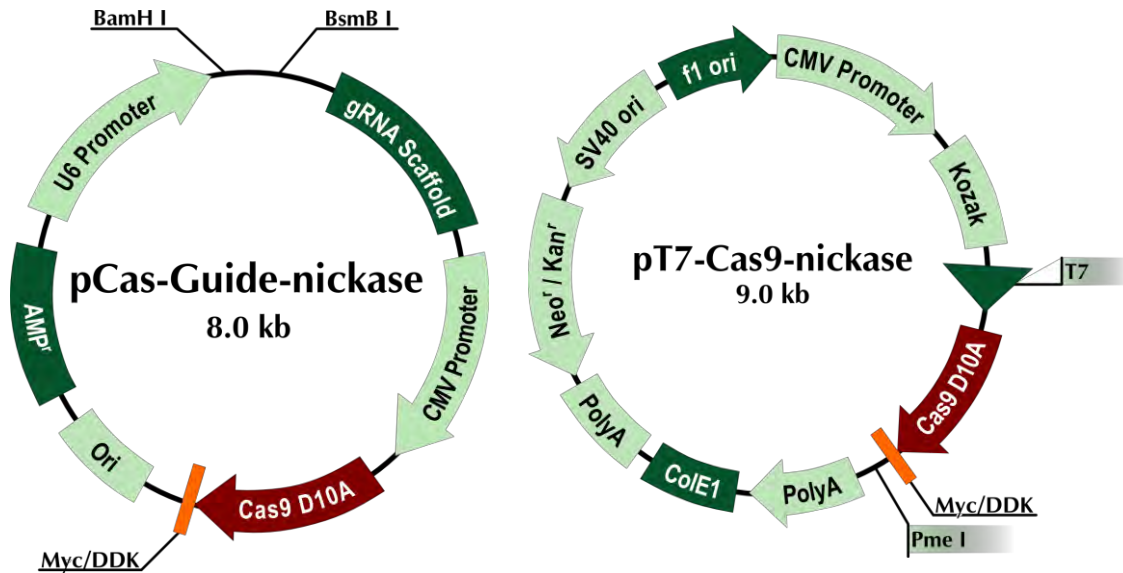

IX. CRISPR scramble controls, pCas-Scramble (SKU GE100003) and pCas-Scramble-EF1A-GFP (SKU GE100021)


**Package contents**

- One (1) vial of circular plasmid DNA, pCas-Scramble (SKU: GE100003) or pCas-Scramble-EF1a-GFP (SKU GE100021), 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis
- Application Guide

**Product description**

Using CRISPR technology for genome editing, you need a negative scrambled control, a 20bp scrambled sequence cloned in CRISPR vectors. After transfecting the CRISPR scrambled control into cells, a guide RNA containing the scrambled sequence (which does not target any sequence) will be produced.

Figure 13. Plasmid maps of pCas-Scramble and pCas-Scramble-EF1a-GFP

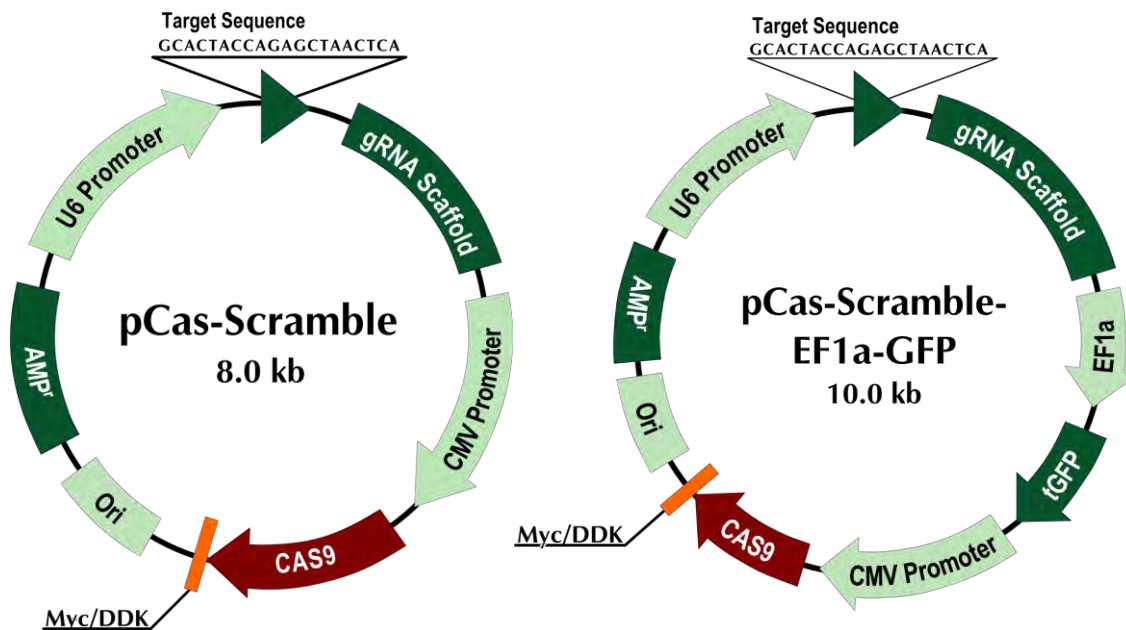


## X. Transgene Knockin via CRISPR at AAVS1 and ROSA26 Loci

### Introduction

Many research projects require inserting a full-length cDNA expression cassette into a host genome for the purpose of stably expressing the encoded protein. In such studies, the integration site of the transgene is very critical. An inserted gene should have the least negative effect to the host cells and the transgene should be genetically stable. Transgene expression via random integration into the genome is subjected to position effects and silencing. In addition random gene insertion might interrupt or activate the neighboring genes. Genomic safe harbor sites are transcriptionally active and transgene insertion does not have adverse effect on the host cell genome. Adeno-associated virus (AAV) integration site (AAVS1) on human chromosome 19 has been accepted as a high gene expression and safe genomic location; for Mouse cells, ROSA26 is proved to be a genomic safe harbor.


OriGene has the largest full length cDNA collection, Human and Mouse. Development of a high efficiency transgene integration kit can facilitate the use of the cDNA collection. Taking advantage of recently discovered CRISPR technology, we developed AAVS1 and ROSA26 safe harbor knockin systems for transgene integration. The AAVS1 or ROSA26 targeting sequence is cloned in the CRISPR all-in-one vector, pCas-Guide; this gRNA/CRISPR vector will generate double stranded genomic break at AAVS1 or ROSA26 locus. The donor vector contains the AAVS1 or ROSA26 homologous arms with the gene of interest in between will be integrated at AAVS1 or ROSA26 site via homologous recombination.

Figure 14. Diagram of transgene insertion at AAVS1 or ROSA26 via CRISPR

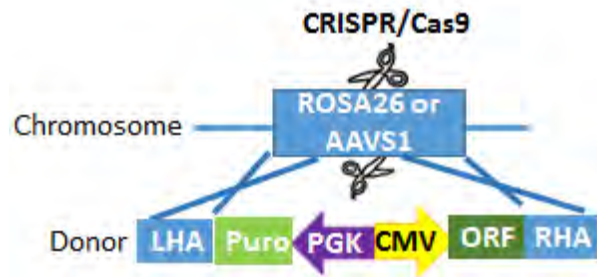

### AAVS1 Transgene knockin via CRISPR

Figure 15. Diagram of how AAVS1 targeted insertion via CRISPR works.

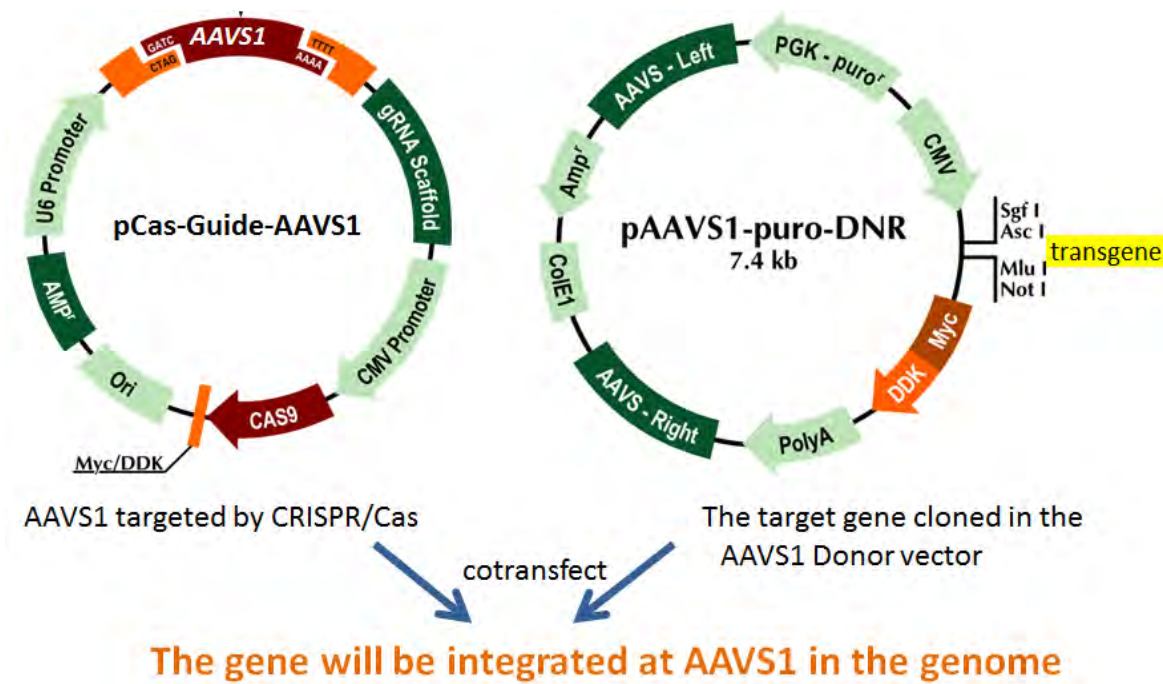

### AAVS1 CRISPR/gRNA vector, pCas-Guide-AAVS1 (SKU GE100023)

#### Package contents

- One (1) vial of pCas-Guide-AAVS1 (SKU: GE100023) plasmid DNA 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit )

The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

### Related Optional Reagents

- Competent E. coli cells
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents


### Related OriGene Products

- pAAVS1-puro-DNR (SKU GE100024)
- pAAVS1-BSD-DNR (SKU GE100035)
- pAAVS1-RFP-DNR (SKU GE100026)
- pAAVS1-Cas9-Puro-DNR (SKU GE100037)
- pAAVS1-Cas9-BSD-DNR (SKU GE100039)
- AAVS1 Transgene knockin vector kit (puro) (SKU GE100027)

### Product Description

pCas-Guide-AAVS1: a plasmid DNA with Cas9 expression and gRNA targeting AAVS1 site. This all-in-one vector will generate a double strand break in human host cells at AAVS1 locus.

Figure 16. Vector map of pCas-Guide-AAVS1.

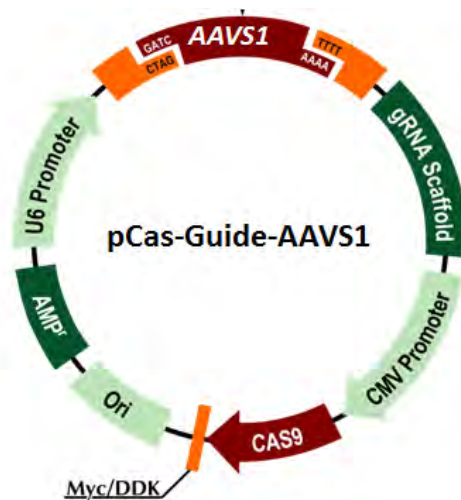



### AAVS1 donor vectors

SKU	Vector Name	Promoter	Cell Selection
GE100024	pAAVS1-Puro-DNR	CMV	Puromycin
GE100035	pAAVS1-BSD-DNR	CMV	Blasticidin
GE100046	pAAVS1-EF1a-Puro-DNR	EF1a	Puromycin
GE100048	pAAVS1-EF1a-BSD-DNR	EF1a	Blasticidin

### Package contents

- One (1) vial of corresponding donor vector, plasmid DNA 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- One (1) Vial of dried 5' (VP1.5) primer (100 picomoles), reconstitute in 10 µL dH<sub>2</sub>O to make a 10 µM solution.
- One (1) Vial of dried 3' (XL39) primer (100 picomoles), reconstitute in 10 µL dH<sub>2</sub>O to make a 10 µM solution.
- Certificate of Analysis

*\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit)*

*The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.*

### Related Optional Reagents

- Competent E. coli cells
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

### Related OriGene Products


- pCas-Guide-AAVS1 (SKU GE100023)
- pCas-Guide-scramble (SKU GE100003)
- pAAVS1-RFP-DNR (SKU GE100026)
- AAVS1 Transgene knockin vector kit (puro) (SKU GE100027)
- AAVS1 transgene knockin vector kit (BSD) (SKU GE100036)
- AAVS1 Transgene knockin vector kit (EF1a-puro) (SKU GE100046)
- AAVS1 transgene knockin vector kit (EF1a-BSD) (SKU GE100048)

## Product Description

AAVS1 donor vectors contain multiple cloning site to clone your gene of interest. The plasmid has AAVS1 left homologous arm and right homologous arm for homologous repair to integrate the transgene cassette into AAVS1 locus. Flanked by the two arms, there is a CMV or EF1a promoter driven expression cassette for a transgene expression and a PGK driven puromycin or blasticidin resistant gene for mammalian selection (Fig. 16). A multiple cloning site downstream of the CMV/EF1a promoter is designed to be compatible to OriGene's precision shuttle vector system so that the ORF clones can be easily shuttled to this integration vector by simply "cut and ligate".

**Note: The AAVS1 donor vectors don't contain your gene of interest; the gene of interest needs to be cloned.**

Figure 17. Vector maps of AAVS1 donor vectors

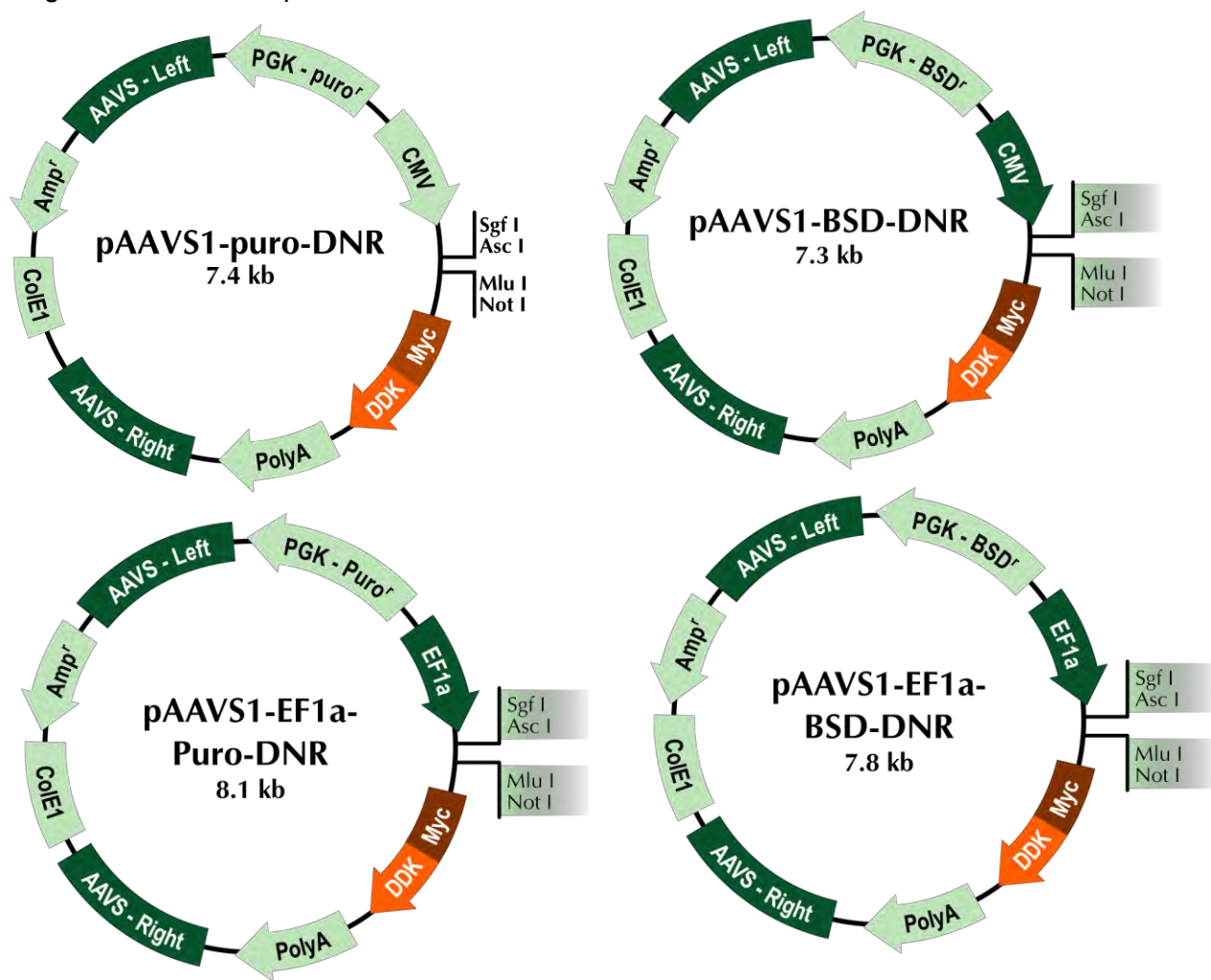

MCS of pAAVS1-Puro-DNR

Kozac  
Consensus

EcoR I Sal I Kpn I RBS Sgf I Asc I

CTATAGGGCGGCCGGGAATTCGTCGACTGGATCCGGTACCGAGGAGATCTGCCGCCGCGATCGCCGGCGCCAGATCT

Mlu I Not I Myc.Tag

CAAGCTTAAGCTAGCTAGCGGACCG ACG CGT ACG CGG CCG CTC GAG CAG AAA CTC ATC TCA GAA GAG  
T R T R P L E Q K L I S E E

DDK.Tag Fse I

GAT CTG GCA GCA AAT GAT ATC CTG GAT TAC AAG GAT GAC GAC GAT AAG GTT TAA ACGGCCGGCC  
D L A A N D I L D Y K D D D D K V Stop

MCS of pAAVS1-BSD-DNR

Kozac  
Consensus

EcoR I Sal I Kpn I RBS Sgf I

CTATAGGGCGGCCGGGAATTCGTCGACTGGATCCGGTACCGAGGAGATCTGCCGCCGCGATCGCC

Asc I Rsr II Mlu I Not I

GGCGCGCCAGATCTCAAGCTTAAGCTTAAGCTAGCTAGCGGACCGACGCGTACGCGGCCGCTC GAG CAG AAA  
E Q K

Myc-DDK tag Pme I Fse I

--- --- GAC GAT AAG GTT TAA ACGGCCGGCC  
D D K V stop

MCS of pAAVS1-EF1a-Puro-DNR & pAAVS1-EF1a-BSD-DNR

Kozac  
Consensus

EcoR I Sal I Kpn I RBS Sgf I

GTCGTGAGCGGCCGGGAATTCGTCGACTGGATCCGGTACCGAGGAGATCTGCCGCCGCGATCGCC

Asc I Mlu I Not I

CGGCGCGCCAGATCTCAAGCTTAAGCTTAAGCTAGCTAGCGGACCGACGCGTACGCGGCCGCTC GAG CAG  
E Q

Myc-DDK tag Pme I

AAA --- --- GAC GAT AAG GTT TAA ACGGCCGGCC  
K D D K V stop

## Cloning a transgene to AAVS1 donor vectors

To clone an ORF insert to the AAVS1 donor vectors, the first step is to select which pair of enzymes. There are a few rare restriction enzymatic sites designed in the multiple cloning site region and c-terminal Myc-DDK tag. Four pairs of restriction enzymes can be used for cloning: SgfI/MluI, Asc/MluI, SgfI/NotI and Asc/NotI. Among them, SgfI/MluI is suitable for over 95% ORF insert. If the ORF inserts are from OriGene's TrueORF clones, they can be easily shuttled from pCMV6-Entry vector to the AAVS1 donor vectors using the above enzyme pairs; the insert will be in frame with c-terminal Myc-DDK tags. In case of none of the enzyme pair can be used, a single enzymatic site can be used for cloning. In such case, screening an insert with the right orientation is needed. If you want to express a native protein without a tag, a termination codon must be introduced at the 3' end of the ORF insert which can be achieved by PCR or site-direct mutagenesis. The following protocol is for shuttling an ORF insert to the donor vectors using SgfI/MluI enzymatic pair.

Digest an ORF insert from TrueORF clone in pCMV6-Entry vector (cat# RC2xxxxxx)

<u>Component</u>	<u>Volume</u>
10X restriction buffer	2 $\mu$ l
Sgf I (10 U/ $\mu$ l)	0.6 $\mu$ l
Mlu I (10 U/ $\mu$ l)	0.6 $\mu$ l
nuclease-free water	13.8 $\mu$ l
<u>TrueORF clone (500 ng)</u>	<u>3 <math>\mu</math>l</u>
Total volume	20 $\mu$ l

Incubate at 37°C for 3 hrs.

2. Digest AAVS1 or ROSA26 donor vectors:

<u>Component</u>	<u>Volume</u>
10X restriction buffer	2 $\mu$ l
Sgf I (10 U/ $\mu$ l)	0.6 $\mu$ l
Mlu I (10 U/ $\mu$ l)	0.6 $\mu$ l
nuclease-free water	14.8 $\mu$ l
<u>AAVS1 or ROSA26 donor vector (200ng)</u>	<u>2 <math>\mu</math>l</u>
Total volume	20 $\mu$ l

\* For the 4% of the clones that have internal Sgf I or Mlu I sites, please use the appropriate combination of restriction sites as recommended by OriGene.

Incubate at 37°C for 3 hrs. Add 0.5  $\mu$ l antarctic phosphatase (units used according to the manufacturer's protocol) to the digestion, and continue to incubate at 37°C for an additional 30 minutes.

3. Purify the digestion using a commercial PCR purification column and elute in 20  $\mu$ l 10 mM Tris.

4. Set up a ligation reaction:

Component	Volume
10 x T4 DNA ligation buffer	1 $\mu$ l
T4 DNA Ligase (4U/ $\mu$ l)	0.75 $\mu$ l
nuclease-free water	3.25 $\mu$ l
digested ORF insert (step 1)	2 $\mu$ l
digested vector (Step 2)	3 $\mu$ l
Total volume	10 $\mu$ l

Incubate the ligation reaction at room temperature for 1 hour.

5. Transform the ligation reaction using high-efficiency competent *E. coli* cells ( $\geq 1 \times 10^8$  CFU/ $\mu$ g DNA) following the appropriate transformation protocol. Plate the transformants on LB-agar plates supplemented with 100  $\mu$ g/ml ampicillin.

6. Pick at least four colonies for subsequent DNA purification and screening. Amplify and purify the selected clone(s) by growing overnight in liquid LB containing the corresponding antibiotics (100  $\mu$ g/ml ampicillin), then isolating the DNA using standard plasmid purification procedures.

Confirm the insert by restriction digestion and/or vector primer sequencing using the provided V1.5 for 5' end sequencing and XL39 for 3' end.

**AAVS1 Transgene knockin vector kit, puro (SKU GE100027, GE100047) or BSD (GE100036, GE100049)**

### Package Contents

1. One (1) vial pCas-Guide-AAVS1 (SKU GE100023), targeting AAVS1 site, 10  $\mu$ g, lyophilized. Reconstitute in 100  $\mu$ l dH<sub>2</sub>O to make a final concentration of 100 ng/  $\mu$ L\*.
2. One (1) vial containing pCas-Guide-scramble (SKU GE100003) negative control, 10  $\mu$ g, lyophilized. Reconstitute in 100  $\mu$ l dH<sub>2</sub>O to make a final concentration of 100 ng/  $\mu$ L\*.
3. One (1) vial containing one of the AAVS1 donor vectors, 10  $\mu$ g, lyophilized. Reconstitute in 100  $\mu$ l dH<sub>2</sub>O to make a final concentration of 100 ng/  $\mu$ L\*.
4. Forward (VP1.5) and reverse (XL39) sequencing primers, 100pmols each, dried onto the bottom of screw cap tubes. Reconstitute each in 10  $\mu$ l dH<sub>2</sub>O to make a 10  $\mu$ M solution. *Primers are for donor vectors only, not for pCas-Guide-AAVS1 or pCas-Guide-scramble.*
5. Certificate of Analysis

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit) .

The cDNA clone is shipped at room temperature, but should be kept at -20°C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

### Related Optional Reagents

- Nuclease free water
- Sgfl and Mlu I
- T4 DNA ligase and buffer
- Competent E. coli cells
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

### Related OriGene Products

- pAAVS1-RFP-DNR (SKU GE100026)
- Genome-wide ORF clones:
- Transfection reagent
- CRISPR/Cas9 products
- Reagents and supplies for immunoblots: user preferred. OriGene has a selection of antibodies and detection reagents.
- DNA purification reagents

### Notice to purchaser

This product is for research use only. Use in and/or for diagnostics and therapeutics is strictly prohibited. By opening and using the product, the purchaser agrees to the following: The plasmids may not be distributed, resold, modified for resale or used to manufacture commercial products without prior written approval from OriGene Technologies, Inc. If you do not agree to the above conditions, please return the UNOPENED product to OriGene Technologies, Inc. within ten (10) days of receipt for a full refund.

### Product description

The AAVS1 transgene knockin vector kits are complete vector kits to knock in your gene of interest in AAVS1 locus in human genome for robust and stable expression.

The kit contains the following three plasmids:

1. pCas-Guide-AAVS1: a plasmid DNA with Cas9 expression and gRNA targeting AAVS1 site (Fig. 15). This two-in-one vector will generate a double strand breakage in the human host cells at AAVS1 locus.
2. pCas-Guide-Scrambled control: the negative gRNA control containing a non-specific gRNA sequence in pCas-Guide vector.
3. AAVS1 donor vector: An AAVS1 repair donor vector containing multiple cloning site to clone your gene of interest. The plasmid has AAVS1 left homologous arm and right

homologous arm for homologous repair to integrate the transgene cassette into AAVS1 locus. Flanked by the two arms, there is a CMV or EF1a promoter driven expression cassette for a transgene expression and a PGK driven puromycin or blasticidin resistant gene for transgene integration selection (Fig. 16). A multiple cloning site downstream of the promoter is designed to be compatible to Origene's precision shuttle vector system so that the ORF clones can be easily shuttled to this integration vector by simply "cut and ligate".

*Note: AAVS1 donor vectors don't contain your gene of interest, it needs to be cloned.*

### Experimental Protocol

1. Clone your gene into the AAVS1 donor vector.
2. Follow the protocols on page: 49-52 in this manual, where pAAVS1-RFP-DNR is used. In your experiment, you can replace your pAAVS1-GOI-DNR with pAAVS1-RFP-DNR.

### pAAVS1-RFP-DNR (SKU GE100026), positive control

#### Package contents

- One (1) vial of pAAVS1-RFP-DNR (SKU: GE100026) plasmid DNA 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis

*\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit)*

*The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.*

#### Related Optional Reagents

- Competent E. coli cells
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

#### Related OriGene Products

- pCas-Guide-AAVS1 (SKU GE100023)
- pCas-Guide-scramble (SKU GE100003)
- pAAVS1-puro-DNR (SKU GE100024)
- AAVS1 Transgene knockin vector kit (puro) (SKU GE100027)
- AAVS1 Transgene knockin vector kit (BSD) (SKU GE100036)
- pAAVS1-Cas9-Puro-DNR (SKU GE100037)
- pAAVS1-Cas9-BSD-DNR (SKU GE100039)

**Product Description**

pAAVS1-RFP-DNR is a positive donor vector with RFP cloned in pAAVS1-puro-DNR vector. This vector can be used as a positive control for CRISPR knockin system when combined with pCas-Guide-AAVS1. pCas-Guide-AAVS1 generate double strand break at AAVS1 site; then pAAVS1-RFP-DNR provides repair template for homologous recombination. The result is RFP-puro expression cassette is inserted in AAVS1 locus.

Figure 18. Vector map of pAAVS1-RFP-DNR


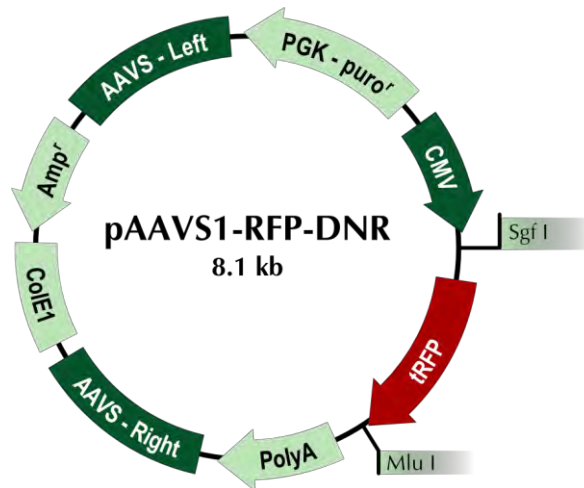


Figure 19. Diagram of CRISPR positive control to knockin RFP at AAVS1 site



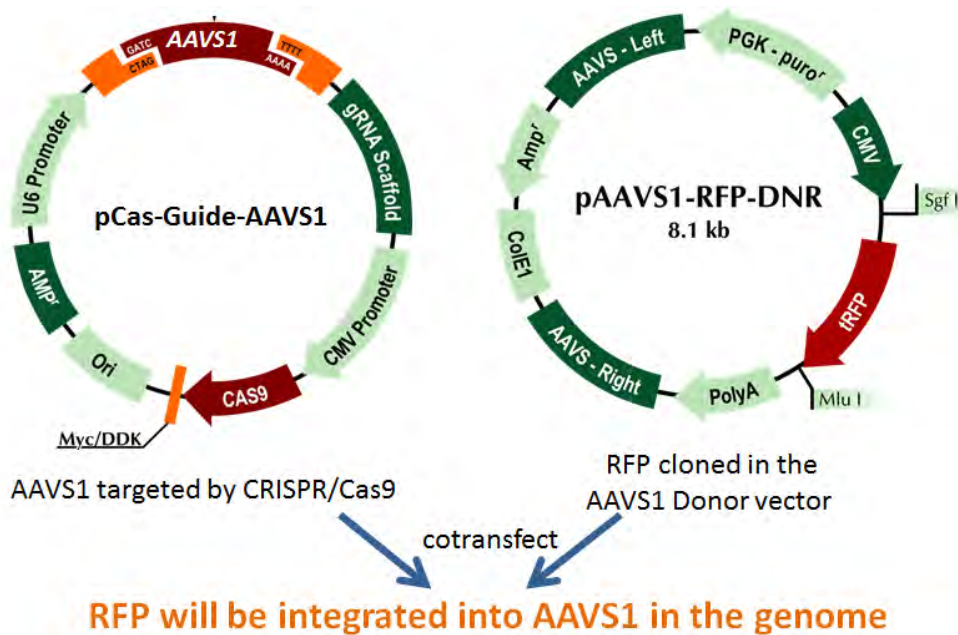


## Experimental Protocol

This protocol is to knockin RFP expression cassette into AAVS1 locus using pAAVS1-RFP-DNR and pCas-Guide-AAVS1; pCas-Guide-scramble is used as a negative control.

A sample protocol listed below is for experiments performed in 6-well plates and using TurboFectin (cat# TF81001) as transfection reagent. If your experiments require a different size of culture plates, just scale up or down the reagents accordingly based on the relative surface area of your plate (Table 1). Different type of cells may need a different transfection reagent; please follow the manufacturer's corresponding protocol. OriGene just launched Virus-like Viromers which are best for difficult-to-transfect cells.

1. Approximately 18-24 hours before transfection, plate  $\sim 3 \times 10^5$  adherent cells in 2 ml culture medium into each well of a 6-well plate or  $\sim 5 \times 10^5$  suspension cells per well to obtain 50-70% confluence on the following day. The amount of cells varies depending on the size of your cells.
2. Transfection in a complete culture medium. Two separate transfections:

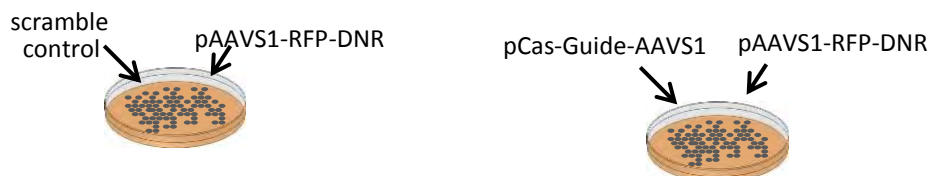

In a small sterile tube, combine the following reagents in the prescribed order. The order of reagent addition is important to achieve the optimal results.

- a. Dilute 1  $\mu\text{g}$  of pCas-Guide-AAVS1 (or scramble control) in 250  $\mu\text{L}$  of Opti-MEM I (Life Technologies), vortex gently. Then add 1  $\mu\text{g}$  of the RFP-donor DNA into the same 250  $\mu\text{L}$  of Opti-MEM I. Vortex gently.
- b. Add 6  $\mu\text{L}$  of Turbofectin 8.0 to the diluted DNA (not the reverse order) and pipette gently to mix completely.
- c. Incubate the mixture 15 minutes at room temperature.

**Note:** We recommend starting with the ratios of Turbofectin 8.0 and DNA listed in table 1; however, subsequent optimization may further increase the transfection efficiency.

**Note:** For your specific cell lines, please choose a transfection reagent that gives you high transfection efficiency and follow the manufacture's protocol.

3. Add the mixture above drop-wise to the cells plated in step 1 (no need to change the media). Gently rock the plate back-and-forth and side-to-side to distribute the complex evenly.
4. Incubate the cells in a 5%  $\text{CO}_2$  incubator.
5. Passage cells around 3 weeks before puro selection or RFP sorting.  
48hr post transfection, split cells 1:10, grow additional 3 days; then split the cells again 1:10. Split cells 7 times in total. Since puromycin resistant gene in the donor vector

contains PGK promoter, RFP is driven by CMV promoter, the plasmid donor DNA before genomic integration will also provide puromycin resistance and express RFP. The reason to passage cells for around 3 weeks before puromycin selection is to dilute out cells containing the donor as episomal form.

### Timelines of genome editing

- ✓ CRISPR targeted gene knockout / knockin--- 1 week post transfection
- ✓ Episomal donor vector dilution with cell passaging--- 3 weeks post transfection

**Note 1.** *Since stable cell selection takes time, you can try to analyze the cells at P2 to detect genomic integration using genomic PCR (Fig 2). When designing primers for genomic PCR, one primer should be outside of the homologous arm region in donor DNA and one primer is in the functional cassette. Please see details in Fig. 20. The amplified PCR fragment is around 1kb. There could be some difficulties doing genomic PCR at this step before selection due to the percentage of edited cells and difficulties of genomic PCR.*

6. Apply puromycin selection or RFP sorting. Since after around 3-week cell passaging, episomal donor DNA is most in most cells, you can use RFP to do cell sorting to enrich edited cells. Another way is to use puromycin selection. Split P5 or P7 cells 1:10, then grow cells directly in the puromycin containing complete media in 10 cm dishes (apply puro selection at P6 might work, but you might get more random donor integrated clones than P8). The dose of puromycin needs to be determined with a kill curve to find out the lowest dose that kills the non-transfected cells completely 4-7 days post selection; the range of puromycin is 1µg/ml to 10µg/ml). Change the media every 2-3 days.

**Note:** *We recommend you still keep growing or freeze some of the transfected cells without selection; just in case, you need to perform the puromycin selection again.*

7. The puromycin resistant cells are ready to be analyzed for genome editing.
  - Use microscope to observe RFP expression.
  - WB with anti-tRFP antibody (cat# [TA150061](#)) to detect RFP expression
  - Genomic PCR to verify the integration of the functional cassette, then directly sequence the amplified genomic fragment using the PCR primers to verify the integration in the genome.

**Note:** *Scramble control and donor vector will also give you some puromycin resistant cells as donor vector alone can randomly integrate into the genome too; however the efficiency should be a lot lower than with a specific gRNA. Therefore, you should get more colonies with gene specific gRNA than scramble control if the gene specific gRNA cleaves efficiently.*

8. Isolate individual cell colonies.
  - Two main methods, limiting dilution and cloning rings / cylinder.
  - 3) Limiting dilution
 

This method is better to be used after puromycin selection. Dilute cells to seed about 1-2 cells/well in 96-well plate, after 1-2 weeks, observe under the microscope and select the wells only containing one cell colony, then further expand them to 6-well plate when they are confluent in the 96-well plate and so on.

#### 4) Cloning rings / cylinder

This method can be used in the same time with puro selection. Seed cells at lower density, such as 5% confluency in a larger cell culture dish, such as 10cm dish, apply puromycin selection when seeding.

Table 1. Recommended starting transfection conditions for Turbofectin 8

Tissue Culture Vessel	Growth area, cm <sup>2</sup> /well	µg of DNA	Ratio of Turbofectin:DNA
96-well plate	0.35	0.1-0.15	3:1
24-well plate	2	0.5-1	3:1
12-well plate	4	1-2.5	3:1
6-well plate	9.5	1-5	3:1
35 mm plate	8	1-5	3:1
60 mm plate	20	2-10	3:1
100 mm plate	60	5-15	3:1

Figure 3. Diagram of knockin process


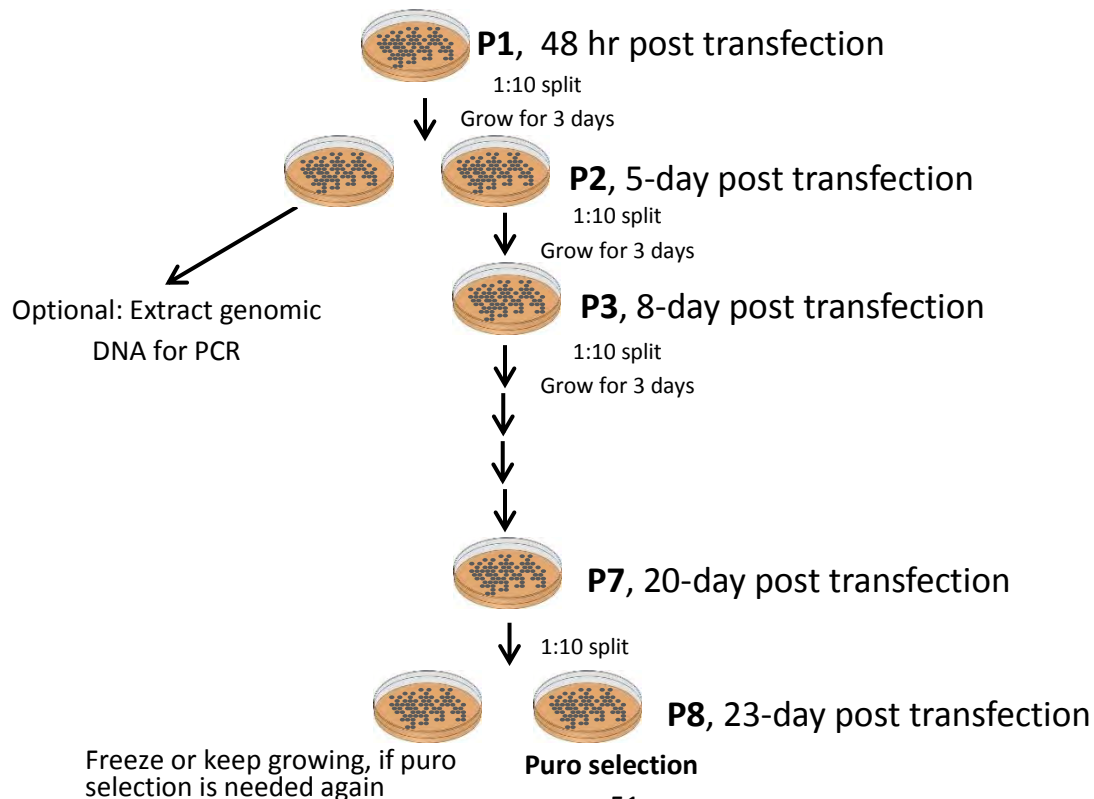


Fig. 4. Diagram of genomic PCR Primer design.

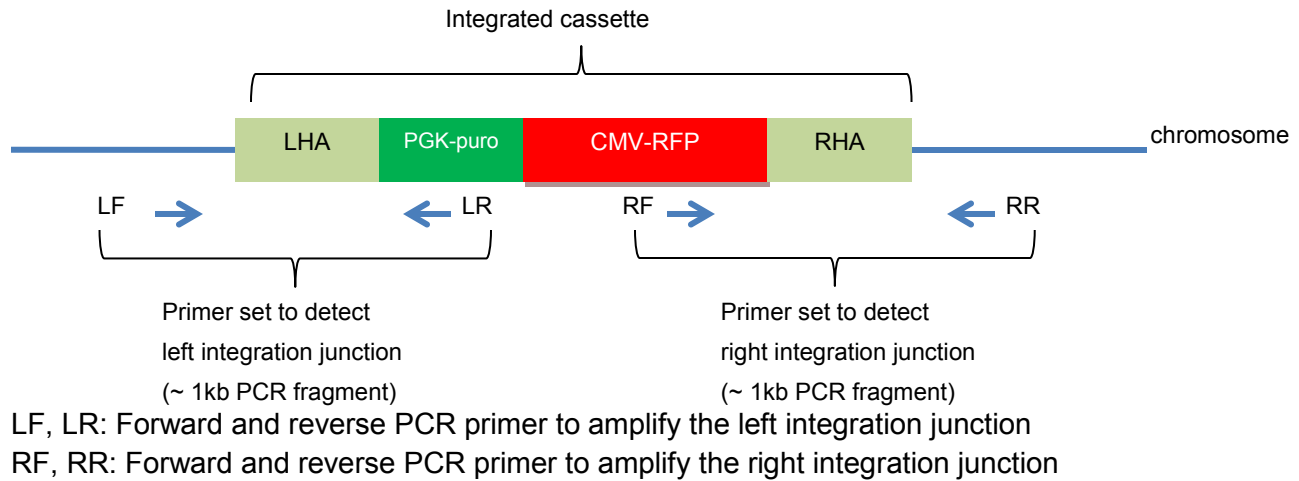

**Cas9 insertion at AAVS1 locus**


**Product description**

For genome editing projects, Cas9 stable cells are needed; to achieve robust and stable Cas9 expression, inserting Cas9 at AAVS1 locus in human cells is preferred. Adeno-associated virus integration site 1 (AAVS1) in human genome is a safe harbor for transgene integration. It is transcriptionally active and transgene expression from this site is robust and stable.

With AAVS1 gRNA targeting construct (GE100023, page 31 of this manual) and AAVS1-Cas9 donor vectors, Cas9 can be inserted at AAVS1 locus. AAVS1 gRNA targeting construct will lead to sequence specific double-strand break at AAVS1 locus in human cells, the AAVS1 homologous sequences in the donor vectors will lead to Cas9-selection cassette being integrated into AAVS1 site vis homologous recombination.

**Cas9 AAVS1 donor vectors**

Figure 20. Plasmid maps of pAAVS1-Cas9-Puro-DNR & pAAVS1-Cas9-BSD-DNR

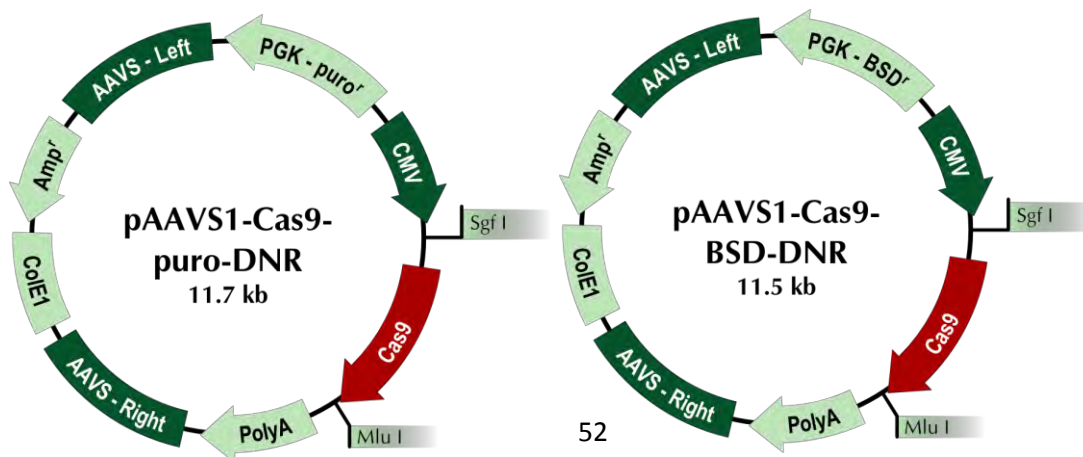


## AAVS1 Cas9 insertion vector kits, Puro (SKU GE100038) and BSD (SKU GE100040)

### Package Contents:

1. One (1) vial of pCas-Guide-AAVS1 (SKU GE100023), 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL\*.
2. One (1) vial of pCas-Guide-scramble (cat# GE100003) negative control, 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL\*.
3. One (1) vial of pAAVS1-Cas9-Puro-DNR (SKU GE100037) or pAAVS1-Cas9-BSD-DNR (SKU GE100039), 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL\*.
4. Certificate of Analysis

Figure 21. Diagram of how Cas9 is inserted at AAVS1 locus in human cells

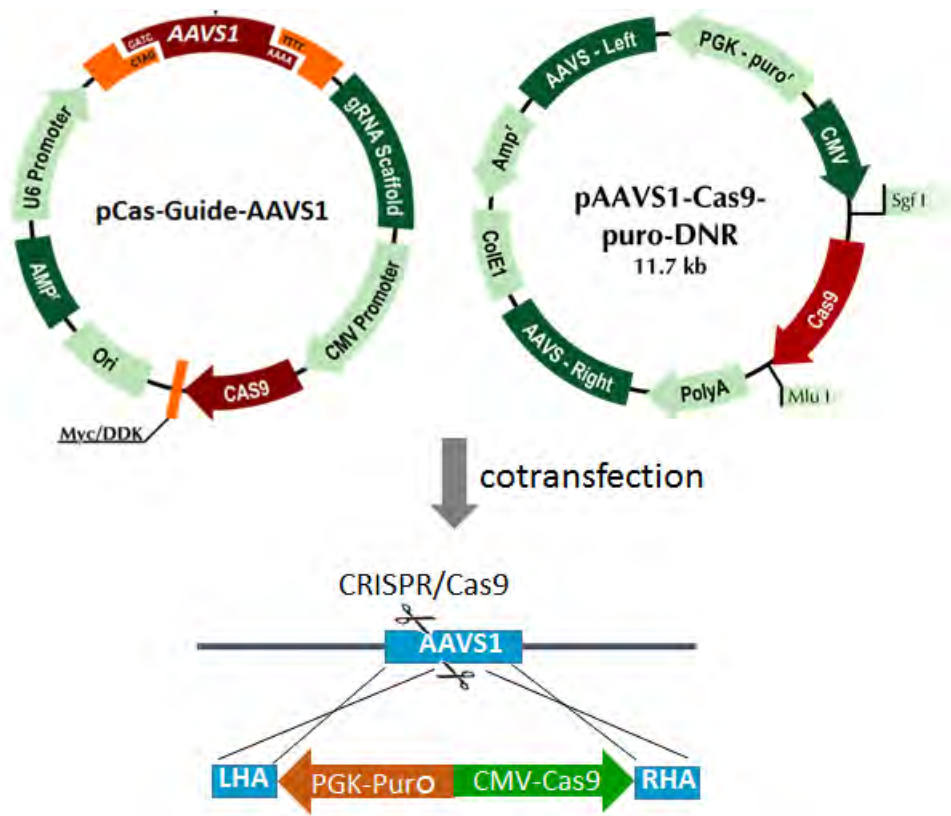


### Experimental Protocol

Inserting Cas9 into AAVS1 site is the same as inserting RFP into AAVS1 site. Please follow the protocols on page 49-52 in this manual, replacing pAAVS1-RFP-DNR with the corresponding Cas9 AAVS1 donor.

## ROSA26 Transgene knockin via CRISPR

Figure 22. Diagram of ROSA26 targeted insertion via CRISPR.

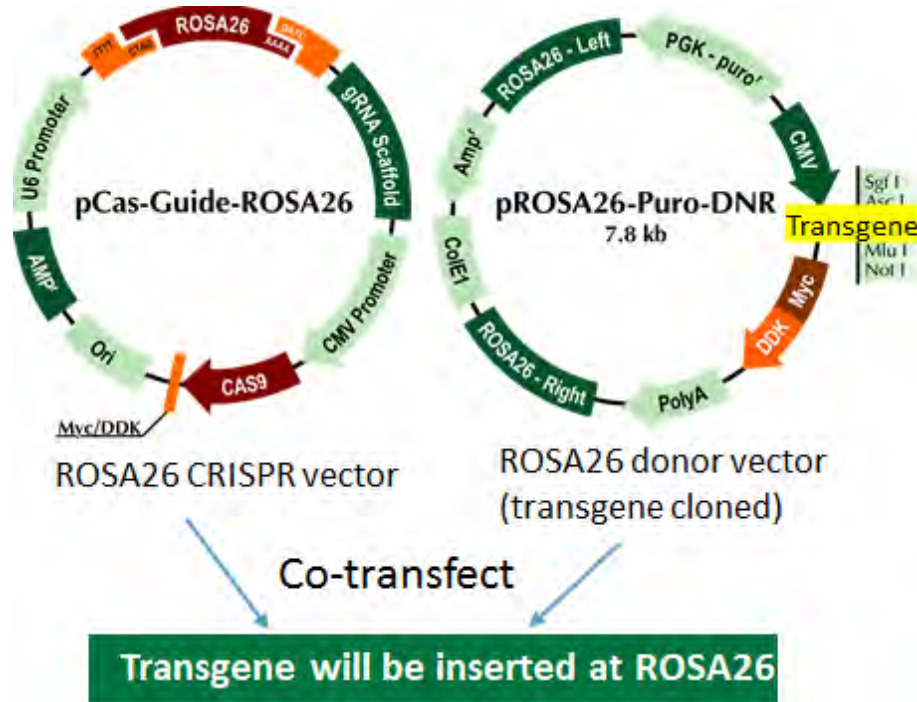

### ROSA26 gRNA/CRISPR construct, pCas-Guide-ROSA26 (SKU GE100050)

#### Package contents

- One (1) vial of pCas-Guide-ROSA26 (SKU: GE100050) plasmid DNA 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- Certificate of Analysis

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit). The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

#### Related Optional Reagents

- Competent E. coli cells
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents


## Related OriGene Products

- pROSA26-Puro-DNR ([GE100051](#))
- ROSA26 transgene knockin vector kit (puro) (SKU [GE100052](#))

## Product Description

pCas-Guide-ROSA26: a plasmid DNA with Cas9 expression and gRNA targeting ROSA26 site. This all-in-one vector will generate a double strand break in human host cells at ROSA26 locus.

Figure 23. Vector map of pCas-Guide-ROSA26.

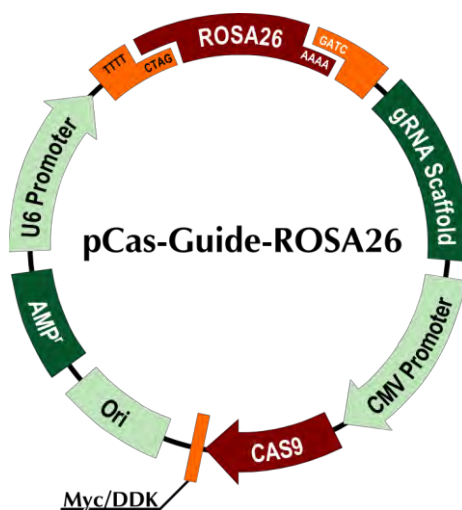

**ROSA26 donor vector, pROSA26-Puro-DNR (SKU GE100051),**

## Package contents

- One (1) vial of pROSA26-Puro-DNR plasmid DNA 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL.
- One (1) Vial of dried 5' (VP1.5) primer (100 picomoles), reconstitute in 10 µL dH<sub>2</sub>O to make a 10 µM solution.
- One (1) Vial of dried 3' (XL39) primer (100 picomoles), reconstitute in 10 µL dH<sub>2</sub>O to make a 10 µM solution.
- Certificate of Analysis

*\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations ([PowerPrep® HP Midiprep Kit](#)). The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.*

## Related Optional Reagents

- [Competent E. coli cells](#)
- LB agar plates with ampicillin, 100 µg/mL

- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

### Related OriGene Products


- pCas-Guide-ROSA26 (GE100050)
- pCas-Guide-AAVS1 (SKU GE100023)
- pCas-Guide-scramble (SKU GE100003)
- pAAVS1-RFP-DNR (SKU GE100026)
- AAVS1 Transgene knockin vector kit (puro) (SKU GE100027)
- AAVS1 transgene knockin vector kit (BSD) (SKU GE100036)

### Product Description

ROSA26 donor vector (pROSA26-Puro-DNR): ROSA26 repair donor vector containing multiple cloning site to clone your gene of interest. The plasmid contains ROSA26 left homologous arm and right homologous arm for homologous repair to integrate the transgene cassette into ROSA26 locus. Flanked by the two arms, there is a CMV promoter driven expression cassette for a transgene expression and a PGK driven puromycin resistant gene for mammalian selection (Fig. 24). A multiple cloning site downstream of the CMV promoter is designed to be compatible to Origene's precision shuttle vector system so that the ORF clones can be easily shuttled to this integration vector by simply "cut and paste".


**Note: pROSA26-Puro-DNR is a donor vector without your gene of interest; the gene of interest needs to be cloned.**

Figure 24. Vector map of pROSA26-Puro-DNR

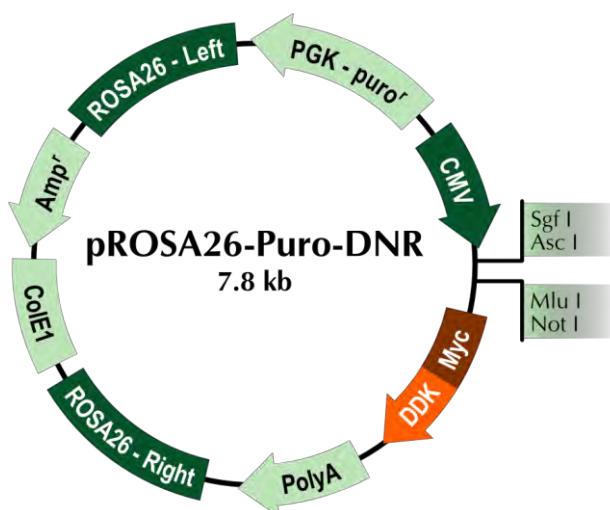



## MCS of pROSA26-Puro-DNR

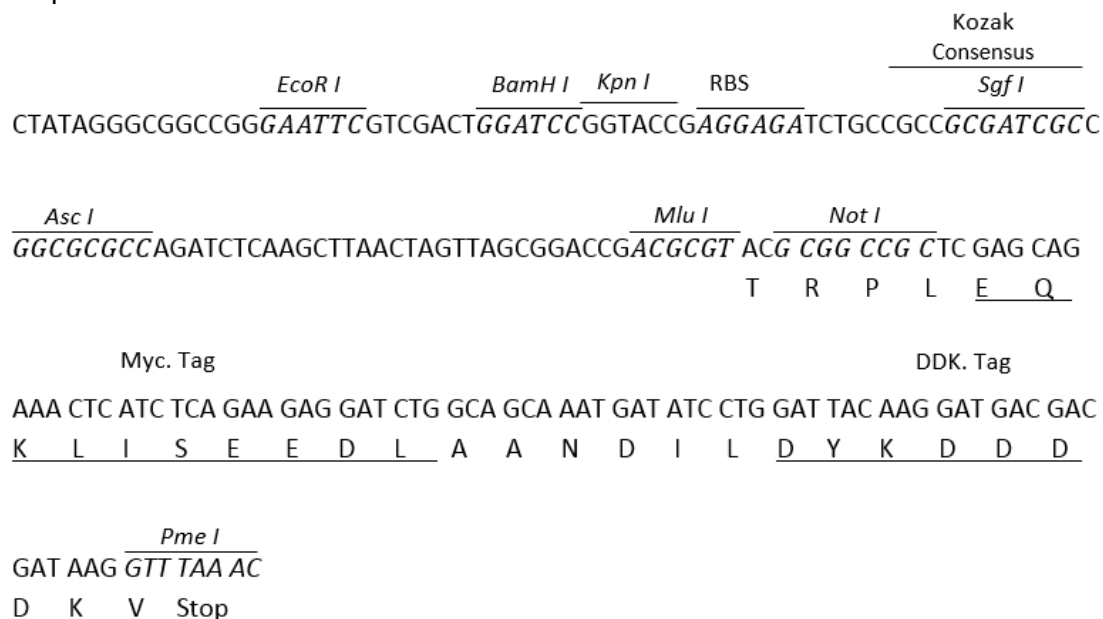

### Cloning a transgene to ROSA26 donor vector

Please follow the protocol on page 44-45 of this manual.

### ROSA26 Transgene knockin vector kit, puro (SKU GE100052)

#### Package Contents

1. One (1) vial pCas-Guide-ROSA26 (SKU GE100050), targeting ROSA26 site, 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL\*.
2. One (1) vial containing pCas-Guide-scramble (SKU GE100003) negative control, 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL\*.
3. One (1) vial containing pRosa-Puro-DNR (SKU GE100051), 10 µg, lyophilized. Reconstitute in 100 µL dH<sub>2</sub>O to make a final concentration of 100 ng/ µL\*.
4. Forward (VP1.5) and reverse (XL39) sequencing primers, 100pmols each, dried onto the bottom of screw cap tubes. Reconstitute each in 10 µL dH<sub>2</sub>O to make a 10 µM solution. *Primers are for donor vector (pROSA26-Puro-DNR) only, not for pCas-Guide-ROSA26 or pCas-Guide-scramble.*
5. Certificate of Analysis

\* OriGene plasmids are purified using ion-exchange columns for high-yield, low endotoxin preparations (PowerPrep® HP Midiprep Kit). The cDNA clone is shipped at room temperature, but should be kept at -20 °C for long-term storage. If properly stored, clones are guaranteed to be stable for 12 months.

### Related Optional Reagents

- Nuclease free water
- Sgfl and Mlu I
- T4 DNA ligase and buffer
- Competent E. coli cells
- LB agar plates with ampicillin, 100 µg/mL
- LB broth (10 g/L Tryptone, 5 g/L Yeast Extract, 10 g/L NaCl. Adjust pH to 7.0 with 1 N NaOH)
- DNA purification reagents

### Related OriGene Products

- Genome-wide ORF clones
- Transfection reagent
- CRISPR/Cas9 products
- Reagents and supplies for immunoblots: user preferred. OriGene has a selection of antibodies and detection reagents

### Notice to purchaser

This product is for research use only. Use in and/or for diagnostics and therapeutics is strictly prohibited. By opening and using the product, the purchaser agrees to the following: The plasmids may not be distributed, resold, modified for resale or used to manufacture commercial products without prior written approval from OriGene Technologies, Inc. If you do not agree to the above conditions, please return the UNOPENED product to OriGene Technologies, Inc. within ten (10) days of receipt for a full refund.

### Product description

The ROSA26 transgene knockin vector kit is a complete vector kit to knockin your gene of interest in ROSA26 locus in mouse genome for robust and stable expression.

The kit contains the following three plasmids:

1. pCas-Guide-ROSA26: a plasmid DNA with Cas9 expression and gRNA targeting ROSA26 site (Fig. 23). This two-in-one vector will generate a double strand breakage in the human host cells at ROSA26 locus.
2. pCas-Guide-Scrambled control: the negative gRNA control containing a non-specific gRNA sequence in pCas-Guide vector.
3. pROSA26-Puro-DNR: A ROSA26 repair donor vector containing multiple cloning site to clone your gene of interest. The plasmid has ROSA26 left homologous arm and right homologous arm for homologous repair to integrate the transgene cassette into ROSA26 locus. Flanked by the two arms, there is a CMV promoter driven expression cassette for a transgene expression and a PGK driven puromycin resistant gene for transgene integration selection (Fig. 24). A multiple cloning site downstream of the CMV promoter is designed to be compatible to OriGene's precision shuttle vector system so that the ORF clones can be easily shuttled to this integration vector by simply "cut and paste".

