Performance guarenteed

# CD8a Monoclonal Antibody (53-6.7), eFluor 450, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Hamster, Fish, Mouse, Human
Host/Isotope	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), eFluor 450, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	53-6.7
Conjugate	eFluor® 450
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage Conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1272198

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.25 µg/test	121 Publications
Functional Assay (FN)	-	2 Publications
Immunocytochemistry (ICC)	-	1 Publication
Immunofluorescence (IF)	-	9 Publications
Immunohistochemistry (IHC)	-	4 Publications
Immunoprecipitation (IP)	-	1 Publication

### **Product Specific Information**

Description: The 53-6.7 monoclonal antibody reacts with the mouse CD8a molecule. CD8a is an approximately 32-34 kDa cell surface receptor expressed either as a heterodimer with the CD8 beta chain (CD8 alpha beta) or as a homodimer (CD8 alpha alpha). A majority of thymocytes and a subpopulation of mature alpha beta TCR T cells express CD8 alpha beta while gamma delta TCR T cells, a subpopulation of intestinal intraepithelial lymphocytes (IELs) and dendritic cells express CD8 alpha alpha. CD8 binds to MHC class I and through its association with protein tyrosine kinase p56lck plays a role in T cell development and activation of mature T cells.

Applications Reported: This 53-6.7 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 53-6.7 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.25  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume


of 100 µL. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Excitation: 405 nm; Emission: 445 nm; Laser: Violet Laser.

Filtration: 0.2 µm post-manufacturing filtered.

## Product Images For CD8a Monoclonal Antibody (53-6.7), eFluor 450, eBioscience™

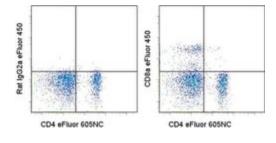

#### CD8a Antibody (48-0081-82) in Flow

Staining of BALB/c splenocytes with Anti-Mouse CD4 eFluor® 605 and 0.125 µg of Rat IgG2a K Isotype Control eFluor® 450 (Product # 48-4321-82) (left) or 0.125 µg of Anti-Mouse CD8a eFluor® 450 (right). Total viable cells were used for analysis.

View more figures on thermofisher.com

### □ 138 References

#### Flow Cytometry (121)

Frontiers in immunology	Species
A New Adjuvant MTOM Mediates Mycobacterium tuberculosis Subunit	Not Applicable
Vaccine to Enhance Th1-Type T Cell Immune Responses and IL-2 <sup>+</sup> T Cells.	Dilution Not Cited
"Published figure using CD8a monoclonal antibody (Product # 48-0081-82) in Flow Cytometry" Authors: Yu Q,Wang X,Fan X	<b>Year</b> 2019
Oncoimmunology	Species
Identification of inhibitors of myeloid-derived suppressor cells activity through phenotypic chemical screening.	Mouse Not Applicable
"Published figure using CD8a monoclonal antibody (Product # 48-0081-82) in Flow Cytometry"	Dilution Not Cited
Authors: Schröder M,Loos S,Naumann SK,Bachran C,Krötschel M,Umansky V,Helming L,Swee LK	Not Cited
	<b>Year</b> 2019

View more Flow references on thermofisher.com

## Immunofluorescence (9)

Nature communications	Species
Oncofetal gene SALL4 reactivation by hepatitis B virus counteracts miR-	Not Applicable
200c in PD-L1-induced T cell exhaustion.	Dilution
"Published figure using CD8a monoclonal antibody (Product # 48-0081-82) in Immunofluorescence"	Not Cited
Authors: Sun C,Lan P,Han Q,Huang M,Zhang Z,Xu G,Song J,Wang J,Wei H,Zhang J,Sun R,Zhang C,Tian Z	<b>Year</b> 2018
PloS one	<b>Species</b>
MHC-I and PirB Upregulation in the Central and Peripheral Nervous	Not Applicable
System following Sciatic Nerve Injury.	<b>Dilution</b>
"Published figure using CD8a monoclonal antibody (Product # 48-0081-82) in Immunofluorescence"	Not Cited
Authors: Bombeiro AL, Thomé R, Oliveira Nunes SL, Monteiro Moreira B, Verinaud L, Oliveira AL	<b>Year</b> 2017

View more IF references on thermofisher.com

### More applications with references on thermofisher.com

IHC (4) IP (1) ICC (1) FN (2)

For Research Use Only, Not for use in diagnostic procedures, Not for resale without express authorization, Products are varranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sells, as set forth in the Production documentation, specifications and for accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year form date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTES, EXPRESS OR IMPLED, ARE GRANTABILTO, IMPLED WARRANTES OF MERCHANTABILTY, FITNESS FOR ANY PARTICULAP VERODS. EXPRESS OR INVERINGENENT, BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTE OR REPLACEMENT OF OR REPLAY. DATE OF MERCHANTABILTY, FITNESS FOR ANY PARTICULAR URPOSE, COR REPLAY DER NELLY EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING TO REPAIR, REPLACE OR REPLAY DER THE NON-CONFORMING PRODUCTS) ATHER IS NO. DO BILGATION THERE IS NO. DO THER WARRANTES OF REPLAY DER NET ON EXCLUSIVE REMEDY OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unloss otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo thereapeutic uses, or any type of consumption

3