

## Technical Data Sheet

## BV510 Mouse Anti-Human CD38


## Product Information

|                         |   |
|-------------------------|---|
| <b>Material Number:</b> | 563251  |
| <b>Alternate Name:</b>  | T10; ADP-ribosyl cyclase 1; Cyclic ADP-ribose hydrolase 1         |
| <b>Size:</b>            | 50 tests  |
| <b>Vol. per Test:</b>   | 5 µl  |
| <b>Clone:</b>           | HIT2  |
| <b>Isotype:</b>         | Mouse IgG1, κ   |
| <b>Reactivity:</b>      | QC Testing: Human   |
| <b>Workshop:</b>        | III 155   |
| <b>Storage Buffer:</b>  | Aqueous buffered solution containing BSA and ≤0.09% sodium azide. |

## Description

The HIT2 monoclonal antibody specifically binds to CD38. CD38 is a 45 kDa type II single-chain transmembrane glycoprotein present on thymocytes, activated T cells and terminally differentiated B cells (plasma cells). Other reactive cells include monocytes, macrophages, dendritic cells and some epithelial cells. The CD38 antigen acts as an ectoenzyme that catalyzes the synthesis and hydrolysis of a Ca<sup>++</sup> mobilizing agent, cyclic ADP-ribose. This intracellular calcium plays an important role in cell signaling pathways. Reports describe CD38 as participating in adhesion with CD31, immunoregulatory functions involving signal transduction leading to cell growth, apoptosis, and differentiation.

The antibody was conjugated to BD Horizon™ BV510 which is part of the BD Horizon™ Brilliant Violet™ family of dyes. With an Ex Max of 405-nm and Em Max at 510-nm, BD Horizon™ BV510 can be excited by the violet laser and detected in the BD Horizon™ V500 (525/50-nm) filter set. BD Horizon™ BV510 conjugates are useful for the detection of dim markers off the violet laser.

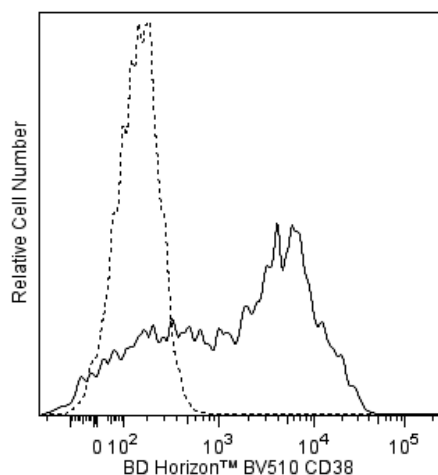

**Flow cytometric analysis of CD38 expression on human peripheral blood lymphocytes.** Human whole blood was stained with BD Horizon™ BV510 Mouse Anti-Human CD38 antibody (Cat. No. 563251; solid line histogram) or with BD Horizon™ BV510 Mouse IgG1, κ Isotype Control (Cat. No. 562946; dashed line histogram). The erythrocytes were lysed with BD Pharm Lyse™ Lysing Buffer (Cat. No. 555899). The fluorescence histograms were derived from events with the forward and side light-scatter characteristics of viable lymphocytes. Flow cytometric analysis was performed using a BD™ LSR II Flow Cytometer System.

## Preparation and Storage

Store undiluted at 4°C and protected from prolonged exposure to light. Do not freeze.

The monoclonal antibody was purified from tissue culture supernatant or ascites by affinity chromatography.

The antibody was conjugated with BD Horizon™ BV510 under optimum conditions, and unconjugated antibody and free BD Horizon™ BV510 were removed.

## Application Notes

## Application

|                |                  |
|----------------|------------------|
| Flow cytometry | Routinely Tested |
|----------------|------------------|

## Suggested Companion Products

| Catalog Number | Name                                | Size   | Clone  |
|----------------|-------------------------------------|--------|--------|
| 554656         | Stain Buffer (FBS)                  | 500 ml | (none) |
| 562946         | BV510 Mouse IgG1, κ Isotype Control | 50 µg  | X40    |
| 555899         | Lysing Buffer                       | 100 ml | (none) |

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

|               |              |               |              |              |                         |
|---------------|--------------|---------------|--------------|--------------|-------------------------|
| United States | Canada       | Europe        | Japan        | Asia Pacific | Latin America/Caribbean |
| 877.232.8995  | 800.979.9408 | 32.53.720.550 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995         |

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.

For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD



## Product Notices

1. This reagent has been pre-diluted for use at the recommended Volume per Test. We typically use  $1 \times 10^6$  cells in a 100- $\mu$ l experimental sample (a test).
2. Source of all serum proteins is from USDA inspected abattoirs located in the United States.
3. An isotype control should be used at the same concentration as the antibody of interest.
4. Please refer to [www.bdbiosciences.com/pharmingen/protocols](http://www.bdbiosciences.com/pharmingen/protocols) for technical protocols.
5. Caution: Sodium azide yields highly toxic hydrazoic acid under acidic conditions. Dilute azide compounds in running water before discarding to avoid accumulation of potentially explosive deposits in plumbing.
6. For fluorochrome spectra and suitable instrument settings, please refer to our Multicolor Flow Cytometry web page at [www.bdbiosciences.com/colors](http://www.bdbiosciences.com/colors).
7. Brilliant Violet™ 510 is a trademark of Sirigen.

## References

McMichael AJ, Beverly PCL, Gilks W, et al, ed. *Leukocyte Typing III: White Cell Differentiation Antigens*. New York: Oxford University Press; 1987. (Clone-specific: Flow cytometry)

Schlossman SF, Boumsell L, Gilks W, et al, ed. *Leukocyte Typing V: White Cell Differentiation Antigens*. New York: Oxford University Press; 1995. (Biology)

## BD Biosciences

[bdbiosciences.com](http://bdbiosciences.com)

|               |              |               |              |              |                         |
|---------------|--------------|---------------|--------------|--------------|-------------------------|
| United States | Canada       | Europe        | Japan        | Asia Pacific | Latin America/Caribbean |
| 877.232.8995  | 800.979.9408 | 32.53.720.550 | 0120.8555.90 | 65.6861.0633 | 55.11.5185.9995         |

For country contact information, visit [bdbiosciences.com/contact](http://bdbiosciences.com/contact)

*Conditions: The information disclosed herein is not to be construed as a recommendation to use the above product in violation of any patents. BD Biosciences will not be held responsible for patent infringement or other violations that may occur with the use of our products. Purchase does not include or carry any right to resell or transfer this product either as a stand-alone product or as a component of another product. Any use of this product other than the permitted use without the express written authorization of Becton, Dickinson and Company is strictly prohibited.*

*For Research Use Only. Not for use in diagnostic or therapeutic procedures. Not for resale.*

Unless otherwise noted, BD, BD Logo and all other trademarks are property of Becton, Dickinson and Company. © 2011 BD

